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Buheung medical – Korea

No.1 healthcare expert
from **Korea**

Thank you very much for buying this product which
has been strictly tested and verified to be safe and reliable.

Buheung

NGUON SONG VIET IMEX TRADING CO., LTD

Ho Chi Minh: 235-237 9A Str., Trung Son, Binh Hung Commune, Binh Chanh Dist, HCMC, Vietnam

Guarantee: 028 7106 5858.

Ha Noi: 2, Lane 382, Pham Van Dong, Co Nhue 2 ward, Bac Tu Liem district, Ha Noi, Vietnam

Guarantee: 024 37 855 966. HOTLINE: 1800 6977

www.buheung.vn

**USER MANUAL
MASSAGE BUHEUNG MK-310 NEW**

Please read all instructions before using


Appropriate Times ©Method to use

- 1 hour after waking up
- 1 hour after meal
- Before sleep
- After a strenuous exercise
- Use 220V-240V outlet
- Hold the handle tight and massager with your speed
- Place the head on the part , and adjust the speed as you like
- For the effective massager,start from the arm ,leg,back, and Side in order
- Drink a glass of water after massage to flush your body properly

Caution

- Do not use the massager more than 20minutes at a time .
- Machine protector is installed ,so when the machine is operated over 20 minutes ,it will stop automatically.
(it is not mean out of order ,switch has to be off and 30minutes later can be re-operated.)
- Do not use the massager without attachments.
- If the supply cord is damaged, it must be replaced by the manufacture or its service agent or a similarly qualified person in order to avoid a hazard.
- The appliance is not intended for use by young children or infirm persons without supervision.
- Young children should be supervised to ensure that they do not play with the appliance.

Speed Control

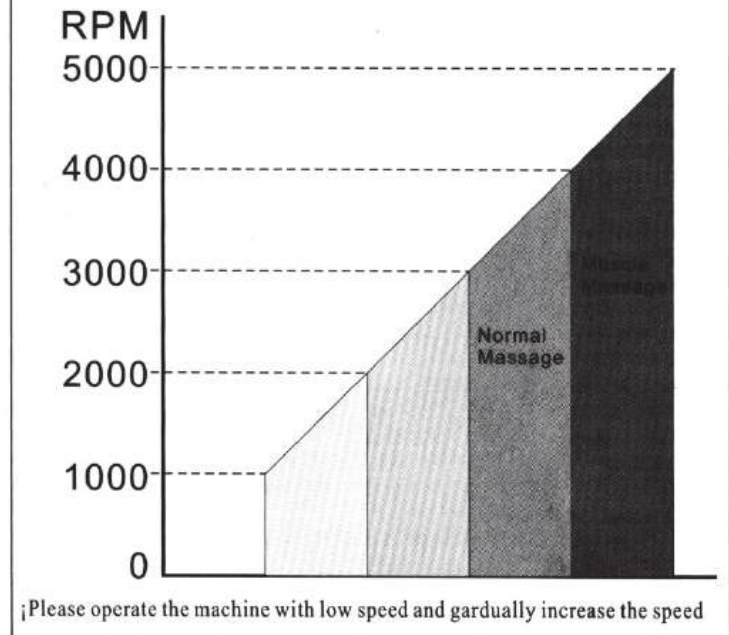Product Specifications

- Volts:220V~240V
- Rated power:25w
- Frequency:50Hz

DANGER

When using a electrical appliance like a Body Massager basic precautions should be followed ,including the following: read all instructions before using the massager

To Reduce the risk of electric shock:

- Always unplug massager from the electrical outlet immediately after using and before cleaning.
- Do not move or store the massager where it can fall or be pulled into a tub or sink
- Do not place in or drop into water or other liquid . Do not reach for a appliance that has fallen into water. Unplug immediately.
- Do not use while bathing or in shower.

WARNING

To Reduce the risk of burns, fire, electrical shock or injury to person:

- Do not operate under a blanket or pillow, or any place where ventilation may be insufficient.
- Close supervision in necessary when used on, by or near children, invalids, or disabled persons.
- Massager should never be left unattended when plugged in .unplug the massager from the outlet when not is use and before cleaning.
- Never operate massager if it has a damaged cord or plug, if it in not working.
- Never drop or insert any object into any openings.
- Use massager only for its intended use as described in this manual. -Do not use attachments not recommended by the manufacture.
- Do not carry massager by the cord or pull the cord sharply.

- Never operate where aerosol products are being used or where oxygen is being administered.
- Do not use outdoors.
- Connect massager el to a properly grounded outlet only.
- To disconnect, turn switch to the off position, then remove from the outlet.
- Keep cord away from heated surface.

PRECAUTIONS

The massager should not be used for:

- Swollen or inflamed areas, or skin eruptions (fractures, bruises & cuts,etc.)
- Do not use on unexplained calf pain, consult your physician.
- Thyroid gland (just below adam's apple:
- Varicose veins
- Eyes, ears, testicles, or ovaries.

- Massage should be comfortable and agreeable,
Whenever in doubt, consult your physician.

- **Pregnant Woman-Special Warning**
Do not use in lower abdomen, lower back, front abdominals, or near ovaries.
Only upper back
Please consult your physician if you have any questions

Please keep the machine like picture when you are not using it

- 1. Please use it under the standard voltage that preset for this product.
- 2. Don’ t use it when the shell was broken lest being injured.
- 3. Don’ t use on the injured or rankle skin .
- 4. Pregnant women and malignant tumor sufferer can’ t use it
- 5. Please don’ t use it in any wet environment.
- 6. Please pull out plug from socket after using it, don’ t pull out the power wire directly or pull out the plug with wet hand.
- 7. Please don’ t allow the children using or playing lest being injured.
- 8. If the power wire was broken or there are some abnormality, please turn off the power(forbid changing by you rself, please sent to our company’ s repair center.
- 9. Please pull out plug from socket when power off, lest it would cause an unforeseeable accident when recov er the power.
- 10.Please turn the switch to ‘OFF’ when finish massage, and pull out plug from socket.
- 11.If you are not a professional repair person, please don’ t open/dismantle this product by yourself.

Trouble shooting

Troubles	Possible reason	Solution
It doesn’ t work after connecting power	The switch is OFF.	Turn on the switch.
	Whether the plug is not correctly insert the socket	Check socket, make sure the plug insert into socket .
The sound become loud suddenly	Whethe pressure is too strong	It will work after reducing pressure
It stop working suddenly	Working for a long time, it stop automatically due to thermal protection..	Cut off power and let it cool down in 30min at least

NEEDLE STICK:
Used to stimulate and tap acupuncture points on the body, specially on the palm, sole of the foot and thighs as you get an acupuncture treatment.

CAUTERY STICK:
Designed for fat removal and used to massage on the belly, waist, thigh and shoulders.

BLOOD LINE MASSAGER STICK:
Used on acupuncture sports and prostate spot .

4 HEAD MASSAGER:
Used on the waist, shoulder, back.

POINT MASSAGER:
Used to massage points of the body.

GEM ARTICULATION STICK:
Used to massage on the joints (articulation) of wrists, ankles, and elbows.

SOFT BRUSH TYPE STICK:
Used on the neck, shoulders, and muscle.

MASSAGER STICK:
Used on weak skin important acupuncture sports, and the thighs.

CUTE MESSAGE STICK:
Used to massage on the joint of body, also used on the waist,shoulde,back and thigh.

FLAT MASSAGER STICK:
Used to massaga muscle,suit for belly, whole back,thigh...etc.