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Thank you very much for buying
this product which has been strictly
tested and verified to be safe and reliable.

Buheung

**PRODUCTS ARE IMPORTED AND DISTRIBUTED EXCLUSIVELY
BY TNHH TM XNK NGUON SONG VIET STORE**

Ho Chi Minh: 235-237 Road No.9A Trung Son residential area,
Binh Hung commune, Binh Chanh district. Ho Chi Minh City.
Hotline : 18008379 | Service Centre : 028 7106 5858

Ha Noi (1st floor): 382, Pham Van Dong , Co Nhue 2 ward,
Bac Tu Liem district, Ha Noi. Hotline : 024 37855966

www.buheung.vn

USER MANUAL

NECK MASSAGER

MK-331

No.1 Health Care from **Korea**

PRODUCT DESCRIPTION

► Product introduction

Smart Neck massages adopts with human body cervical physiological curvature of the popular type circular design, Low-frequency pulse, magnetic effect, hot moxibustion, highfrequency vibration, circular traction, five-in-one synchronous health care to form a composite energy field.


► It can be used in this following cases:

1. Long time exposure before screen.
2. Vervous, insomnia and dreaminess.
3. Stiffness on neck and Muscle fatigue.

► Product Function:

1. Unique 3D grip ergonomics designed.
2. Using the latest 3D intelligent bonding technology, intelligent identification of human skin.
3. Multi-Function, Tntelligent control, operation is simple.
4. There are 6 pulse modes: cupping, tapping, scraping, shiat-su, acupuncture and relaxation.
5. Low-frequency pulse, magnetic effect, hot moxibustion, highfrequency vibration, circular traction, five-in-one synchronous health care to form a composite energy field.
6. Type-C powetr jack.
7. Wireless control.
8. 15 level intensities.
9. 15 minutes timer.
10. Blue light will be flashing with long "DI" every two minutes when out of power. It will be hearing "DI DI DI" warning when pres 'ON/OFFE" button.

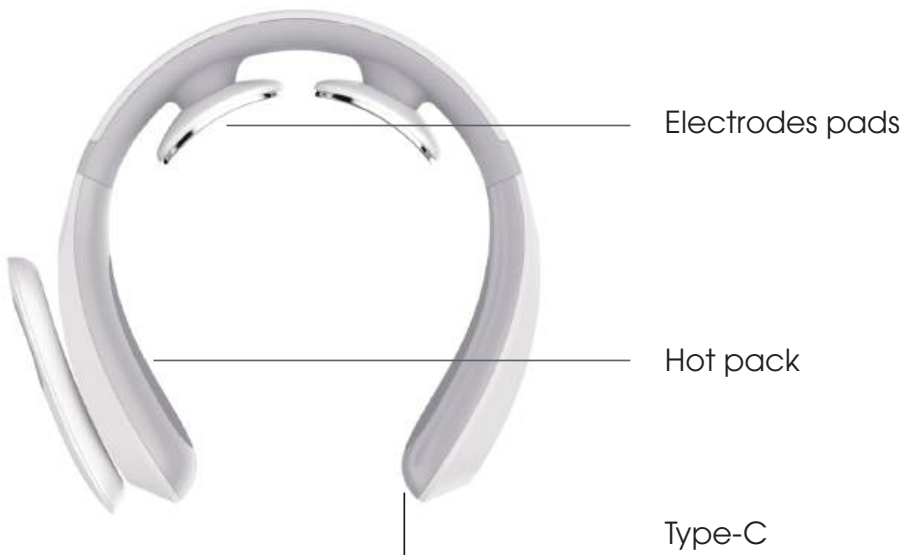Front views


Controller


Stereograms

Device operation

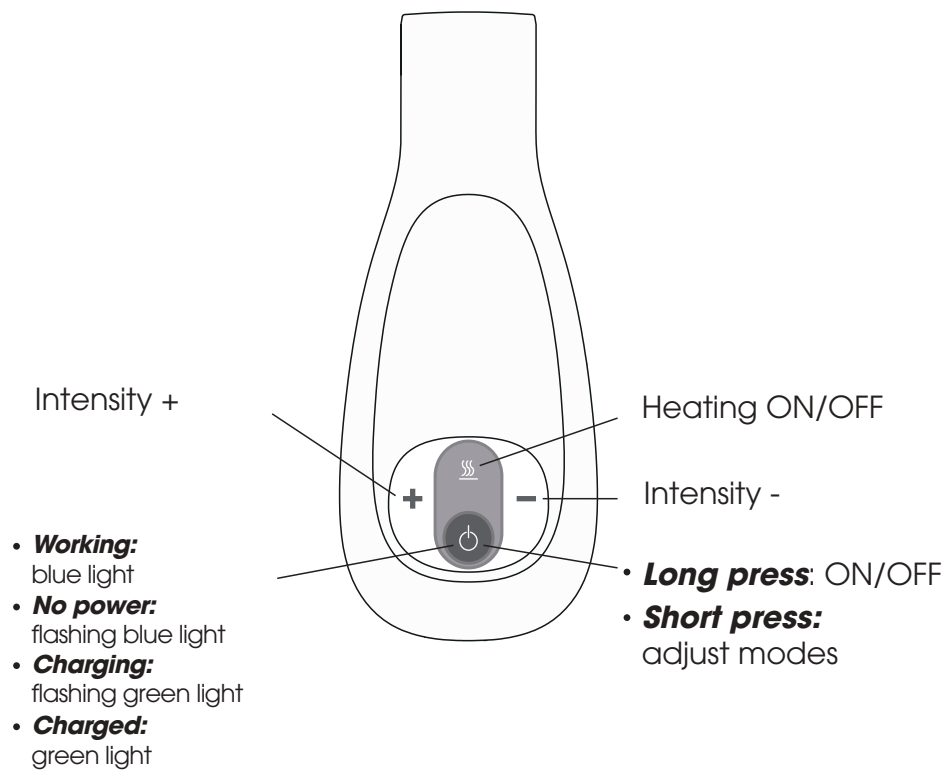Attention

Please keep electrodes pads close to skin.
If not, it will be warning three "DI" when pressing "+"

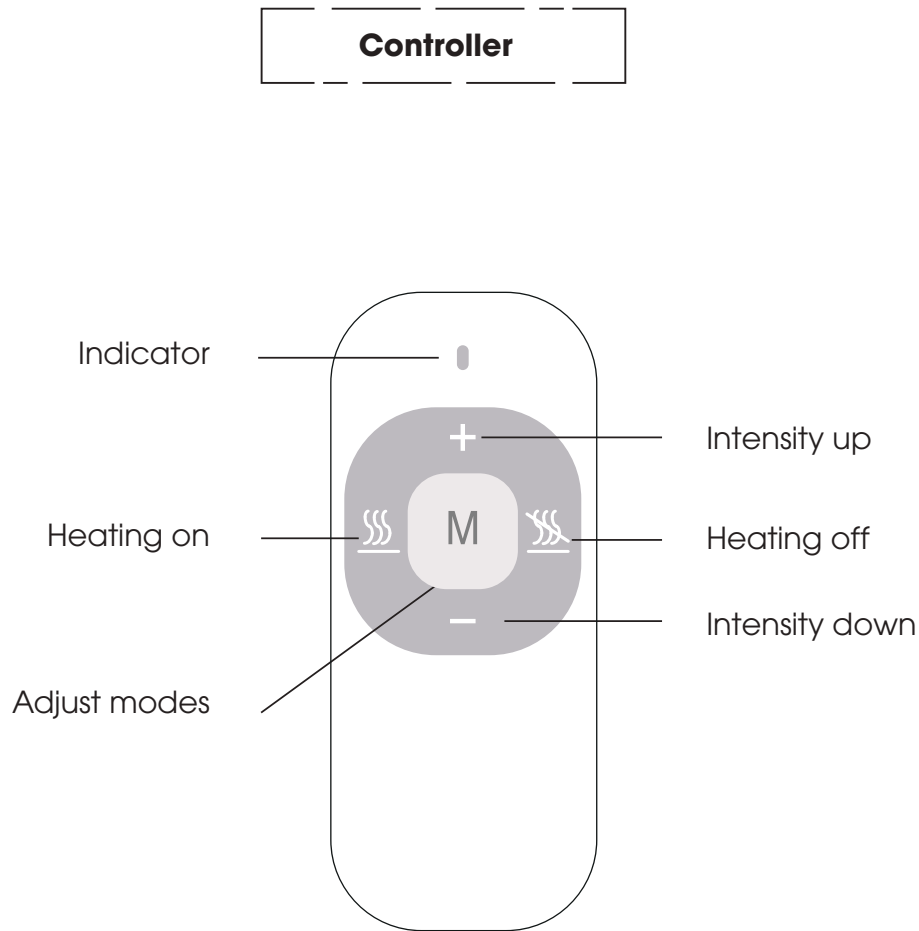Note: Battery CR2032

► Charging:

1. Please charge when out of power, adapter is 5V, 1A. Indicator turns green with flash when charging, it will turn green when full charge. Charging time about 3 hours.
2. Using time about 5 hours when full charging, If massaging with heat function always, using time is 50 minutes.
3. Suggest to charge this device one time per month if do not use for a long time.

► CHÚ Ý:

1. Please do not use this device with other medical machine together, Such as pacemaker, heart lung machine and Medical electronic device, electrocardiograph etc. If patient is using both high frequency surgical equipment and this device, it may hurt skin or damage this device. If using device near the short-wave, microwave treatment equipment (1 meter), EMS function may not be stable.
2. If people who are Pregnant, sensitive skin, heart disease, abnormal blood pressure, malignant tumors, cerebrovascular patients, acute disease, or being treated, please consult your doctor before using this device.
3. Do not use in bath or sleep.
4. Do not use for children.
5. Do not use in patients who are skin sensory impairment or for heat insensitive.
6. Please turn off the device before adjust position or remove device to other place.
7. Please adjust pulse intensity which is able to accept the comfort and appropriate intensity.
8. Please do not disassemble this device, avoiding collision and water.
9. Do not wash this device with water.

TREATMENT

► Treatment:

In order to avoid environmental pollution, please do not throw away this device or accessories when they are aging. Please comply with the local environmental protection requirements.

► FAQ:

- **Q:** What is the cause of some tingling on neck when using electrical pulse function?
- **A:** The reason is that people's skin is relatively dry, electrodes are not close enough to the neck, and cannot stand EMS function. Please keep wet on neck and close enough to skin. Pulse intensity should be low to high. Please change another mode to massage if cannot stand pulse function.
- **PS:** Pulse function is imitating traditional Chinese acupuncture, it is not electric leakage.
- **Q:** Why neck skin in rashes?
- **A:** Neck skin rashes cause skin too dry or not close to neck skin when using pulse function. We suggest to keep neck skin wet and tight close to device, rashes will be gone one or two days.
- **Q:** How long does it take for fully charge? How many times can be use?
- **A:** Built-in Li-ion battery, Charging time is about 3 hours. It can be used in 5 hours only Pulse function, it can be used in 50 minutes if always turn on heating function together.
- **Q:** Why does pinch hurt?
- **A:** Ergonomic design, Perfect fit for different neck sizes. "The product adopts ABS plastics, it can protect neck; Please wear this device correctly and maintain correct posture when using, We suggest to massage two times per day, 15 minutes per time.

► Cleaning and maintenance:

1. Please clean it with soft wet cloth or towel to wipe, cannot use Volatile preparation.
2. Please keep it in a dry and ventilated place.
3. Do not put it in direct sunlight, high temperature and damp place.

► Technical parameters:

Impulse frequency	1-1000Hz
Pulse width	20-400μs
Heat moxibustion temperature	+38°C ~ 48°C

(25°C environment temperature)

► Specifications:

Name	Smart Neck massager
Model	MK-331
Input	DC5V/1A
Rated power	5W
Battery	Built-in Li-ion battery 3.7V/1200MAH
Controller battery	CR2032

► Transport and storage:

Temperature range	-20°C ~ -55°C
Relative temperature range	≤ 80%
Atmospheric pressure range	50kPa ~ 106kPa

(Storage in a dry, ventilated place)

► Element of harmful substances in products:

Name	Hazardous substance					
	Pb	Hg	Cd	Cr(VI)	PBB	PBDE
Device	×	●	●	●	●	●
Accessories	×	●	●	●	●	●
Packing	×	●	●	●	●	●

(This form is made in accordance with SJ/T 11364)