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MASSAGE CHAIR

BUHEUNG **MK-5800**

Scan the Zalo QR code to
activate online warranty

Thank you for choosing and supporting BUHEUNG KOREA!

We hope that you are satisfied with our product and service, and that we have the opportunity to serve and accompany you and your family.

If you have any feedback or reviews regarding the quality of our product and service, please contact our hotline at **1800 8379** so that we can receive and promptly address them.

In addition, please scan the **Zalo OA** code to activate the electronic warranty, receive additional gifts, and enjoy **6 months of warranty** and **24/7** support from BUHEUNG.

Thank you once again and we are happy to have take care of you and your family's health!

Buheung

FEATURE SPECIAL

There are **6 pre-installed** automatic massage programs on the control panel that can be activated with just one touch: **Relaxation, Recovery, Deep Tissue, Flexibility, Pain Reduction, and Sleep.**

Despite being in the affordable price range, the chair's product quality is top-notch, **featuring a 4D massage mechanism with a silicon cover.**

Voice control allows users to activate an automatic massage program with just two commands.

There are 5 intensity levels for the massage that can be adjusted to accommodate most family members.

The airbag system is distributed throughout the chair and covers the feet, calves, buttocks, hips, arms, and shoulders to stimulate blood flow throughout the body.

The chair can be set up in **sleep mode**, a feature exclusive to high-end models, which keeps the chair in a reclined position for resting after a massage session.

CONTENT

- I. Content
- II. Safety notes
 - 1. Operating environment
 - 2. People who should not use this device
- III. Preparing for use
 - 1. Installing cushion
 - 2. Adjust the seat
 - 3. Move the massage chair
 - 4. Position the massage chair
 - 5. Power on
- IV. Names and functions of main components
- V. Introduction to controller functions
- VI. Introduction to massage functions
- VII. Voice control
- VIII. Maintenance
- IX. Specifications

Note: The manufacturer reserves the right to improve the product without prior notice.
The appearance of the product may vary depending on the actual product.

II. SAFETY NOTES

1. Operating environment:

- Do not use or store the device in a damp environment such as a bathroom, sauna, hot spring, etc. to avoid mold, electric shock, fire or damage.
- Do not use the device in a place with a temperature above 40°C.
- Do not use this device outdoors.
- Do not place the device in a place with too much dust, oil, or lack of oxygen.
- Do not place the device near a gas stove, heater or in a place where direct sunlight shines on it, as it will cause the internal parts of the massage chair to overheat and cause fire or explosion.
- Install and use the massage chair on flat surfaces such as floors.
- To avoid electric shock, please confirm if the socket is grounded. Please consult with a technical staff member when installing the grounding wire into the socket.
- Do not smoke while using the massage chair.
- Do not step on the rollers or airbags. Pay attention when entering and exiting the massage chair, and if anything falls out, please reinstall it in its original position.
- When the airbags are being inflated, do not retract your hands or feet, as this will cause the airbags to deflate and return to their original state.

2. People who should not use this device:

The following people should consult a doctor before use (using the device may cause discomfort, lead to accidents and pose a danger to the body):



- Patients with malignant tumors.
- Patients with severe heart disease.
- Patients with severe skin diseases.
- Patients with serious bone and joint diseases.
- Pregnant women.
- Individuals undergoing medical treatment as prescribed by a doctor, especially those experiencing abnormal physical symptoms.
- Individuals with abnormal spinal curvature or spine abnormalities due to accidents or illness.
- Individuals with high fever who need rest and quiet.
- Individuals with abnormal joint symptoms.
- Individuals with medical devices such as "pacemakers".
- Individuals with internal organ diseases (such as stomach, intestinal, dysentery, hepatitis,...).
- Individuals with impaired blood circulation in the legs.
- Heat-emitting devices must be used with caution by individuals who are not sensitive to heat.
- This product is not recommended for individuals with physical or mental disabilities.
- If children are allowed to use this product, an adult should be present to supervise and guide them.
- Children should be supervised to ensure that they do not misuse this product.

Do not allow individuals under 18 years of age or individuals who cannot control their behavior (such as those with mental illness, intellectual disability, physical disabilities, etc.) to use the massage chair alone or touch other parts of the chair. This can cause electric shock, injury to the user, or malfunction of the machine.

III. PREPARING FOR USE

1. Installing cushion

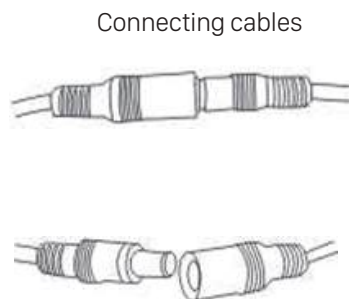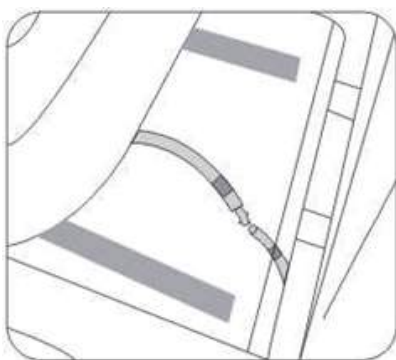Connect/remove the cable

2. Adjust the seat

To adjust the back of the seat, press and hold the manual adjustment switch to adjust the tilt angle of the seat back

Note:


It is necessary to ensure there is enough space in front of and behind the massage chair when adjusting the backrest and footrest.

3. Move the massage chair

Raise both armrests before moving the massage chair


Note:

- Do not have anyone sit on the massage chair while moving it.
- Do not move the chair on a surface that can be easily damaged

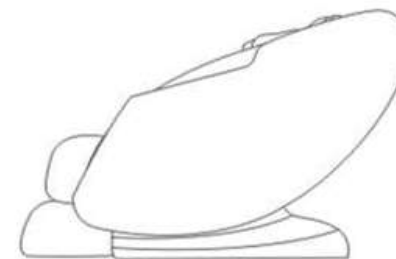4. Position the massage chair


- Please ensure that there is enough space around the massage chair when preparing for a massage.
- The location of the massage chair should have enough space for the back of the chair to recline and the footrest to be raised.
- Make sure there is **at least 80cm** of empty space behind and **at least 50cm** in front of the massage chair.
- Ensure that there are no obstacles in front of or behind the massage chair.
- To avoid damaging the floor, please place a mat underneath the massage chair.

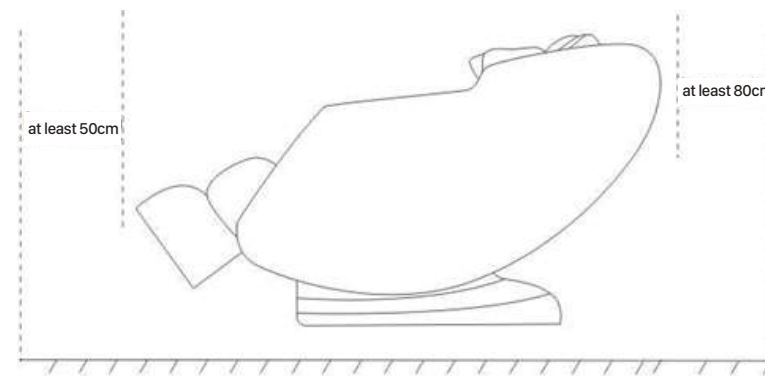
5. Power on

Plug the power cord into an electrical outlet and turn on the power switch.


- ① Please check and confirm if the power cord is being pressed by the massage chair.
- ② Check if the grounding wire is properly connected.
- ③ Before plugging in the power plug, check all buttons of the massage chair, including the power switch.

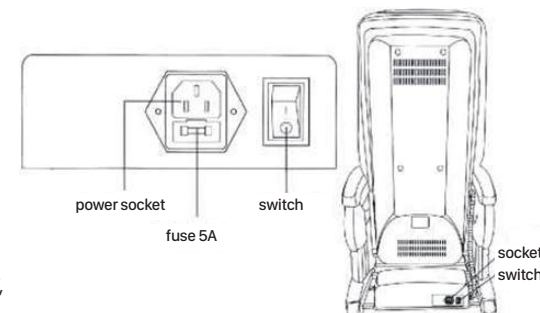IV. NAMES AND FUNCTIONS OF MAIN COMPONENTS

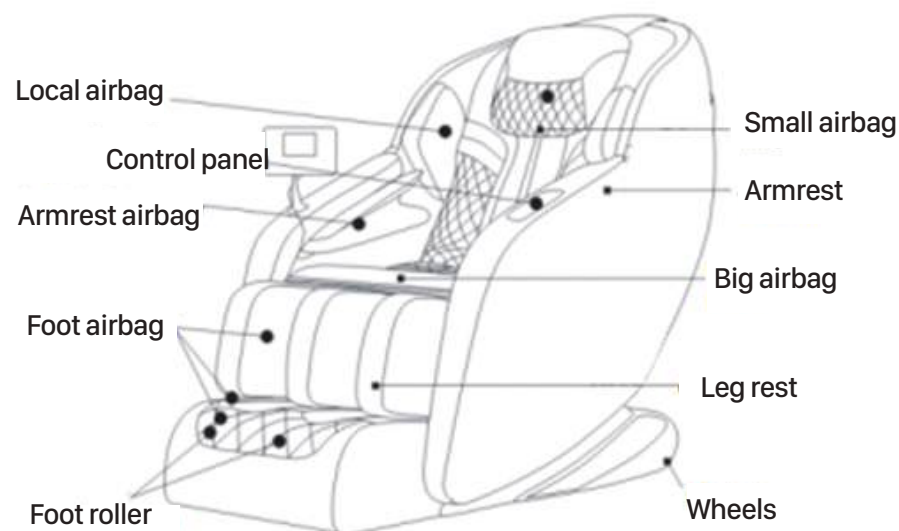V. INTRODUCTION TO CONTROLLER FUNCTIONS

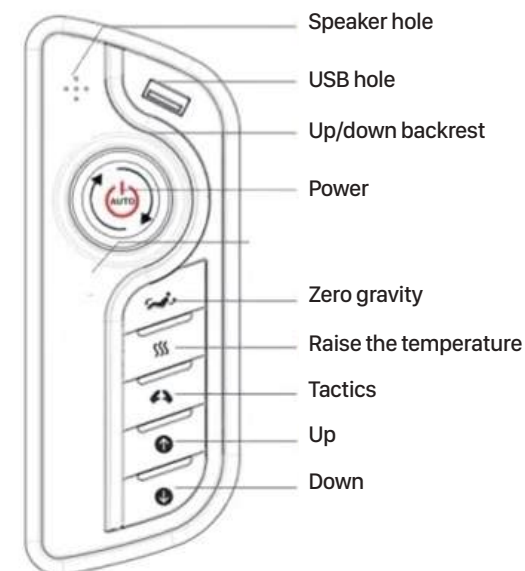VI. INTRODUCTION TO MASSAGE FUNCTIONS

1. Power button:

Press this button  to turn on/off the massage function of the chair. After plugging in the power, press this button to return the chair to its initial position and automatically perform body positioning. After completion, the default first program - "Relaxation Massage" will run and the default running time is 20 minutes.

2. Automatic massage program: Depending on the massage time and physical condition, users can choose one of the five massage modes: "Relaxation", "Recovery", "Flexibility", "Pain Relief", "Sleep" (See LCD screen instructions for more details).

3. Adjustable backrest:

- Press this button  and the backrest will slowly rise while the footrest lowers. Once raised to the desired angle, release the button to lock the backrest and footrest in place.
- Press this button  and the backrest will slowly lower while the footrest raises. Once lowered to the desired angle, release the button to lock the backrest and footrest in place.

4. In the automatic program, press the button **CHUYÊN SÂU** to select massage functions such as "kneading", "punching", "pressing", "shiatsu", "tapping", "acupressure", etc.

5. **Adjusting position:** In manual mode, the massage position for the back can be adjusted up or down by pressing the button **LÊN**

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Voice control only available for Vietnamese

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12. **"Spot" mode:** Press **TỪNG PHẦN** to turn on/off the spot massage mode.

13. **Heat:** Press **NHIỆT** to turn on/off the heat function.

Introduction to shortcut keys:

- Power button: Press and hold for three seconds to turn on/off the machine.
- Zero-gravity button: Press this button to turn on/off the zero-gravity mode.
- Heat button: Press this button to turn on/off the heating mode.
- Up button: Press this button to control the massage of the upper back.
- Down button: Press this button to control the massage of the lower back.

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*Note: Wake up machine with command "Xin chào"

*Note: Voice control only available for Vietnamese

VIII. MAINTENANCE

1. Always keep the massage chair clean and free of dust, and place it in a dry location.
2. If the massage chair is not going to be used for an extended period of time, cover it with a cloth.
3. Do not place the massage chair in direct sunlight or in a location with high temperatures, as this may cause the leather to fade and age quickly.
4. When cleaning the massage chair, please turn off the switch and unplug the power cord.
5. Do not use cleaning agents containing organic solvents such as benzene to clean the massage chair, as this will cause the leather to discolor, age quickly, and develop cracks.
6. If the power cord is damaged, contact the manufacturer or maintenance department to have it replaced.

IX. SPECIFICATIONS

Product name: Massage chair Buheung MK-5800

Model: MK-5800

Rated voltage: 220V~

Rated frequency: 50HZ

Rated power: 120W

Rated operation time: 20 minutes