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No.1 healthcare expert
from **Korea**

Thank you very much for buying this product which has been strictly tested and verified to be safe and reliable.

Buheung

NGUON SONG VIET IMEX TRADING CO., LTD

Ho Chi Minh: 235-237 9A Str., Trung Son, Binh Hung Commune, Binh Chanh Dist, HCMC, Vietnam
Guarantee: 028 7106 5858.

Ha Noi: 2, Lane 382, Pham Van Dong, Co Nhue 2 ward, Bac Tu Liem district, Ha Noi, Vietnam
Guarantee: 024 37 855 966. HOTLINE: 1800 6977

www.buheung.vn

USER MANUAL
MASSAGE CHAIR BUHEUNG MK-5300

Please read this manual carefully before you use it

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:
Read all instructions before using (this appliance).

DANGER-To reduce the risk of electric shock:

- Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- Never use pins or other metallic fasteners with this appliance.
- Carefully examine the covering before each use. Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.
- **WARNING**--To reduce the risks of burns, fire, electric shock or injury to persons:
 - An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
 - Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
 - Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 - Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
 - Do not carry this appliance by supply cord or use cord as a handle.
 - Keep the cord away from heated surfaces.
 - Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
 - Never drop or insert any object into any opening.
 - Do not use outdoors.
 - Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
 - To disconnect, turn all controls to the off position, then remove plug from outlet.
 - Keep children away from extended foot support.
 - Connect this appliance to a properly grounded outlet only. See Grounding Instructions.
 - Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
 - Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.


SAVE THESE INSTRUCTIONS

GROUNDING INSTRUCTIONS

This product must be grounded. If it should malfunction or break down, grounding provides a path of least resistance for the electric current to reduce the risk of electric shock.

This product is equipped with a cord with an equipment-grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances. Danger-Improper connection of the equipment-grounding conductor can result in a risk of electric shock.

Check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product - if it will not fit the outlet, have a proper outlet installed by a qualified electrician.



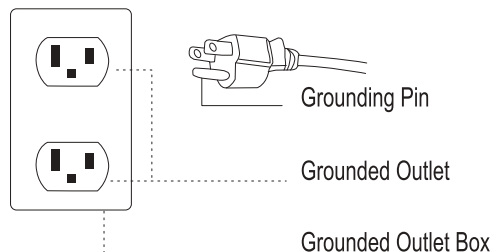
IMPORTANT SAFETY INSTRUCTIONS

OPERATING SURROUNDINGS

● OPERATING ENVIRONMENT

- Do not use or store the massage chair under damp conditions:
 - i.e. bathroom, sauna, spa, to avoid any mold development, electric shock, fire or mechanical failure.
- Don't use the massage chair outdoors
- Don't use the massage chair under dusty and temperature varying conditions.
- Don't place the massage chair near heater(s), stove or direct sunlight to avoid inner components overheating which may lead to electrical fires and shock.
- Make sure to place the massage chair on flat level ground surface
- The chair should be placed on carpet, wood or hard surface flooring
- Make sure the power plug is properly grounded to avoid shock and injury.
- Don't use the machine while smoking

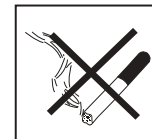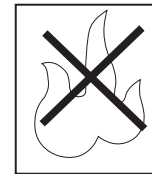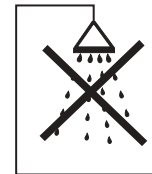FOLLOWING CONDITIONS REQUIRING CONSULT FROM A MEDICAL DOCTORS BEFORE USE

- Suffered from malignant tumor
- Serious heart disease
- Serious skin disease
- Bone disease
- While pregnant or menstruating
- On rest by the doctor
- Spinal condition including injuries associated to the back and spine
- Hypertension and cramping
- Arthritic bone disease
- Any implanted medical device such as pacemaker
- Suffering from internal conditions (i.e. gastritis, enteritis, diarrhea, Zepatitis)
- Serious blood circulation trouble
- Consciousness issues such as vertigo and other equilibrium matters

WARNING

Children under the age of 18 years old or with physical or mental condition are not permitted to use this product unattended by an able adult and with the permission of a medical doctor. These precautions are taken to prevent injury and damage to the product

IMPORTANT SAFETY INSTRUCTIONS

For safety reasons and in order to clearly indicate potential dangers should the product be used improperly, instructions have been marked in the following manners and should be strictly noted

⊘ CAUTIONS

Non-authorized technicians are not permitted to dismantle or repair this product to avoid fire, electrical shock or serious injury

Do not clean or wipe down this product with anything wet or moist when this symbol is shown.

When storing or during maintenance, be sure to disconnect the power plug from socket. Do not detach the plug with wet hands

Make sure to turn the power switch to the "off" position and detach the power plug when cleaning or storing the product

**During cleaning or when moving the chair, make sure to detach power cord from the wall and chair to avoid fire or electrical shock
Do not place the machine on a unstable, uneven and unlevelled surface to avoid damage, insulation failure or electrical shock**

When storing the machine, make sure to disconnect power plug from the outlet and chair to avoid fire caused by electrical surges.

Clean power plug periodically to avoid buildup leading possible fire

Do not use it at any unsafe area or for any other purposes than personal use. Not for commercial use.

Do not store the chair in areas that are not climate controlled. Must be in room temperature with very low humidity.

Keep the power cord unobstructed at all time, i.e. not twisted, knotted, over bent or pulled and do not place anything heavy on it

IMPORTANT SAFETY INSTRUCTIONS

⊘ CAUTIONS

Do not make any modifications or change components by an unauthorized person

Do not dismantle the machine by unauthorized person to avoid electrical shock, short-circuit or fire

Do not store it near a fire, heater or flammable substances It is one of the major reasons of causing product damage.

Should any abnormality occur(i.e. burning smell), stop using immediately and detach power plug from outlet.

please contact your local seller or authorized service agent for maintenance to avoid electrical fire and injury.

Only proper rated voltage must be used . Stop using immediately if power cord or power plug is damaged or not tightly fitted into the outlet or the base of the chair

**Be sure to keep children and pets away while operating to avoid any injury
Do not place the chair in direct sunlight at all times**

If the supply cord is damaged,it must be replaced by the manufacturer or an authorized technician

FCC Radiation Exposure Statement:

This equipment complies with FCC radiation exposure requirement.
 The device can be used in portable exposure condition without RF striction.

Changes or modifications not approved by the party responsible for compliance could void the users authority to operate the equipment.


This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.


If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/ TV technician for help.


FCC Radiation Exposure Statement:

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. In order to avoid the possibility of exceeding the FCC radio frequency exposure limits, Human proximity to the antenna shall not be less than 20cm (8 inches) during normal operation.

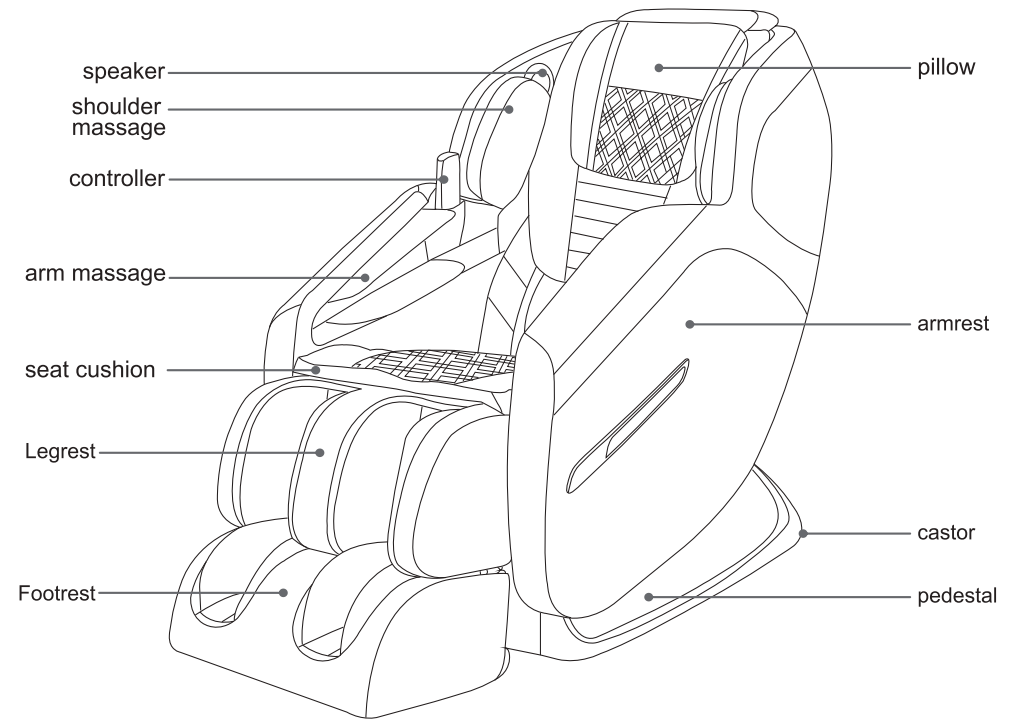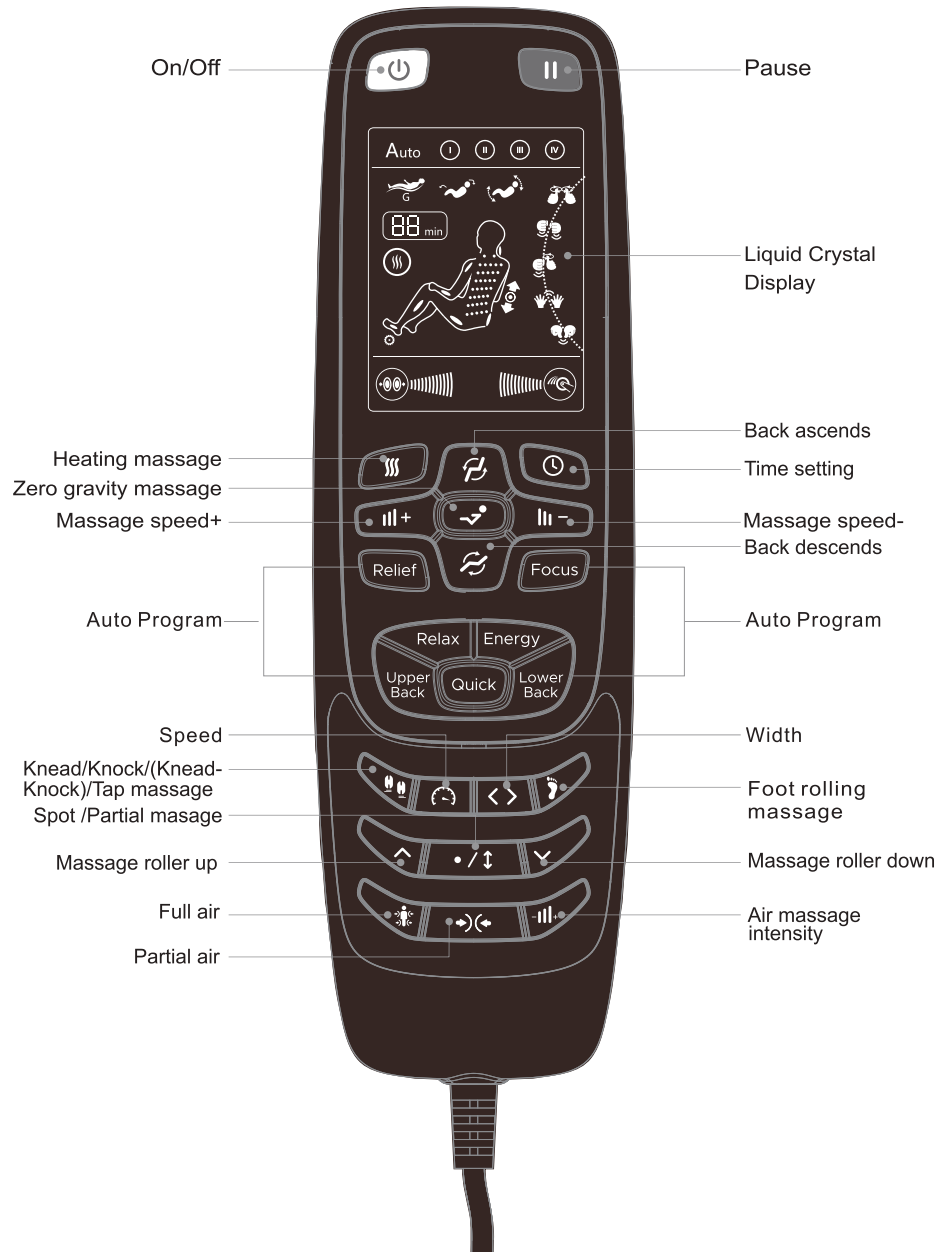Power up the chair

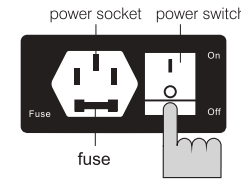Plug in the power cable to the back of the chair and then connect the power cable to the power supply. Press the power switch to the “I” position to power up the chair.

Set time for massage

On/Off

Sit in chair and press the “On” button on the remote controller to switch on the chair. The chair will be in spa mode.


Time setting

Press the "time setting" button to adjust massage time. You can adjust time to maximum 30 minutes. The unit of increment is 5 minutes.

Adjust backrest cushion and leg rest cushion

Set the backrest and leg rest to your desired angle using the “Back ascends” and “Back descends” keys respectively. Press the appropriate button for each adjustment and release your hold once you have set the position .

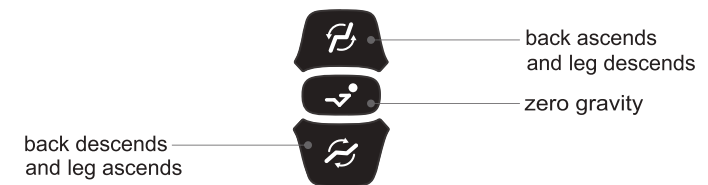Stop or end massage

Pause

The default setting for the massage is 15 minutes after which the massage will automatically stop. The chair settings will remain unchanged from your adjustments. Press “Pause” to stop the massage and press it again to continue. Press the “Off” button and the chair will stop and go into standby mode.

Massages – Automatic programs

The automatic massage programs comprise kneading, tapping, Shiatsu and air massage.

BUTTON	DESCRIPTION	FUNCTIONS
	Relief massage program	Kneading massage program Angle can be adjusted during massage.
	Focus massage program	Performs a series of slower, deliberate massage techniques to focus on muscle groups starting at the top of the back and moving down. Angle can be adjusted during massage.
	Relax massage program	Kneading massage, air massage and tapping massage. Angle can be adjusted during massage.
	Upper back massage program	Upper back massage and leg massage. Angle can be adjusted during massage.
	Lower back massage program	Lower back massage and leg massage. Angle can be adjusted during massage.
	Energy massage program	Kneading massage, air massage and tapping massage. Angle can be adjusted during massage.
	Quick massage program	Press for a 6 minute demonstration of all available massage techniques.

- * For the above programs, the default time is 15 minutes and the default air massage strength is low. Press “Air massage strength” to adjust the intensity of the air pressure.
- * The automatic full-body massage programs cannot be combined with other programs.

 **Warning:** Only rated supply 120V must be used to avoid damage to electronic components.

Massages – Manual programs

Partial air massage is a massage program which targets different needs and parts of the body.

BUTTON	DESCRIPTION	FUNCTIONS AND ADJUSTMENTS
	 Kneading	Press button once to start kneading massage. Speed is adjustable.
	 Knocking	Press button twice to start kneading-knocking massage. Speed is adjustable.
	 Knead-knocking	Press button three times to start knocking massage. Speed is adjustable.
	 Tapping	Press button four times to start tapping massage. Speed is adjustable.
	 Shiatsu	Press button five times to start Shiatsu massage.
	Massage speed regulation	Press this button to adjust massage speed.

- * For the above programs, the default setting is spot-specific massage with medium width. Speed and width can be adjusted with the respective buttons.
- * The above manual massage programs cannot be combined with other programs.

Manual control description and functions

BUTTON	DESCRIPTION	FUNCTIONS
	Speed adjustment	Three levels: slow, medium and fast.
	Width adjustment	Three levels: narrow, medium and wide.
	Full course/ Partial massage	Press button once to start full course massage. Press button twice to start partial massage. Press button three times to stop.
	Massage roller move up	Press button to move massage roller. Release hold on button to stop.
	Massage roller move down	

BUTTON	DESCRIPTION	FUNCTIONS
	Air pressure intensity adjustment	Air pressure intensity is adjustable. It can only be adjusted in air pressure massage mode. Intensity can be soft, medium or strong.
	Full air	Press once to start air massage for upper body. Press twice to start air massage for lower body. Press three times to start air massage for full body. Press four times to stop.
	Partial air	Press once to start air massage for shoulders. Press twice to start air massage for arms. Press three times to start air massage for leg. Press four times to start air massage for foot. Press five times to stop.
	Rolling massage	Foot rolling control. Press button to start foot rolling massage. Press again to stop.
	Heating function	Waist heating function. Press button to start waist heating massage. Press again to stop.

Speakers function

You can pair up the speakers on the chair with your mobile phone with Bluetooth. The massage chair Bluetooth code is BK8000L. Use you mobile phone Bluetooth setting to search for the massage chair device and auto-pair with it. You can then playback music form your mobile phone through the built-in speakers of the massage chair.

Bluetooth Code: BK8000L.

⚠ Warning: Supervision and guidance are necessary when people with disabilities and mental illness. Please advise with your doctor before using the chair if you have ANY medical conditions. Children under 18 should not use the chair and consult with a doctor for approval if attempted. Do not operate the chair if there are any children or pets near by. Serious injuries may occur if these instructions are not followed.

SAFETY INSTRUCTIONS

1. Should any abnormality occur, shut off and detach power plug from socket.
2. When power cord damaged, it must be changed from special power cord or available from other manufacturer or service department.

⚠ Warning: Do not dismantle, modify or change components by unauthorized person to avoid electric shock or fire.

DAILY CARE AND MAINTENANCE

Storage guide

- ⦿ Do not dismantle any part the chair or make modifications. Serious injury make occur
- ⦿ Do not attempt to repair the machine by unauthorized person to avoid electric shock or fire.
- ⦿ Make sure to shut off and detach the power plug when cleaning to avoid electric shock or hazard.
- ⦿ Keep the massage chair clean and free dust. Do not place the massage chair in a humid area.
- ⦿ If not in use for a long period, please cover it with a cloth to protect from dust.
- ⦿ Do not place the massage chair under sunlight or heating source like a heater or directly by a window that has direct sunlight.

Maintenance

- ⦿ Use a neutral detergent only that is free of acid and alcohol.
- ⦿ Use a neutral detergent to clean the back cushion, head cushion and leg cushion. Let it dry.
- ⦿ With a damp cloth wipe down the plastic parts, pipe and Synthetic leather material of the chair.
- ⦿ Make sure not to use any gasoline, amyl acetate, benzene based cleaners to avoid leather color fading or discoloring.
- ⦿ Do not spray pesticide or aerosols near or on the chair to avoid electric shock or fire.
- ⦿ Please clean with dry cloth to avoid mechanical failure.

Controller maintenance

- ⦿ Do not place the massage chair near heater or fire.
- ⦿ Keep away from humidity to avoid electric shock and damage to the chair.
- ⦿ Please clean dust with soft cloth or small banister brush.
- ⦿ Make sure not to use detergents containing benzene, gasoline, amyl acetate.

TROUBLE SHOOTING AND SPECIFICATION

TROUBLE SHOOTING

Checking is necessary in frequent usage.

Make sure the machine work normally before usage, after not in used for long time.

Abnormality situation	Solution guide
Failure to start and indicator light is off	<ul style="list-style-type: none">⊙ Check the plug is well connected⊙ Check the switch on the chair is turned on in the “I” position
Sound in operation	<p>The following sounds are normal</p> <ul style="list-style-type: none">⊙ The mechanism will produce a sound for kneading or knocking⊙ Sound changes after speed adjustment⊙ Sound from when the rollers go up or down⊙ When the chair does a tapping massage it will make the chair vibrate making a sound⊙ When the air inflates and deflates it will make a swooshing sound⊙ Sound for airbag inflating and venting. Airbags may crackle⊙ Sound for leg adjustment⊙ Sound for motor when the massage chair is loaded

SPECIFICATION

Product name	Massage Chair
Model	Titan TP-Lotus
Power supply	120V~ 60Hz
Rated power	150W
Rated working time	15 min