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**PRODUCTS ARE IMPORTED AND DISTRIBUTED EXCLUSIVELY BY
NGUON SONG VIET GROUP**

Ho Chi Minh : 235-237 Road No.9A Trung Son residential area, Binh Hung commune, Binh Chanh district.
Ho Chi Minh City. **Hotline** : 18006977. **Service Centre**: 028 7106 5858

Ha Noi: 382 Pham Van Dong street, Co Nhue 2 ward, Bac Tu Liem district, Ha Noi, Vietnam
Service Centre: 024 37855966

www.buheung.vn

**MK-5150 MASSAGE CHAIR
PRODUCT MANUAL**

Please read the Manual carefully before installing and using the product

Thank you for using
our company's products and services!

Buheung


Thank you for purchasing your Ultimate Health massage chair
Before operating your new chair, please read this instruction manual carefully
Please pay attention to the important safety warnings and keep this manual for future use.

CONTENTS


CHAIR OVERVIEW
UNPACKING AND INSTALLATION
INSTALLATION
MOVING THE CHAIR
REMOTE CONTROL MANUAL
REMOTE CONTROL KEYS
PRODUCT SPECIFICATIONS

01 CHAIR OVERVIEW

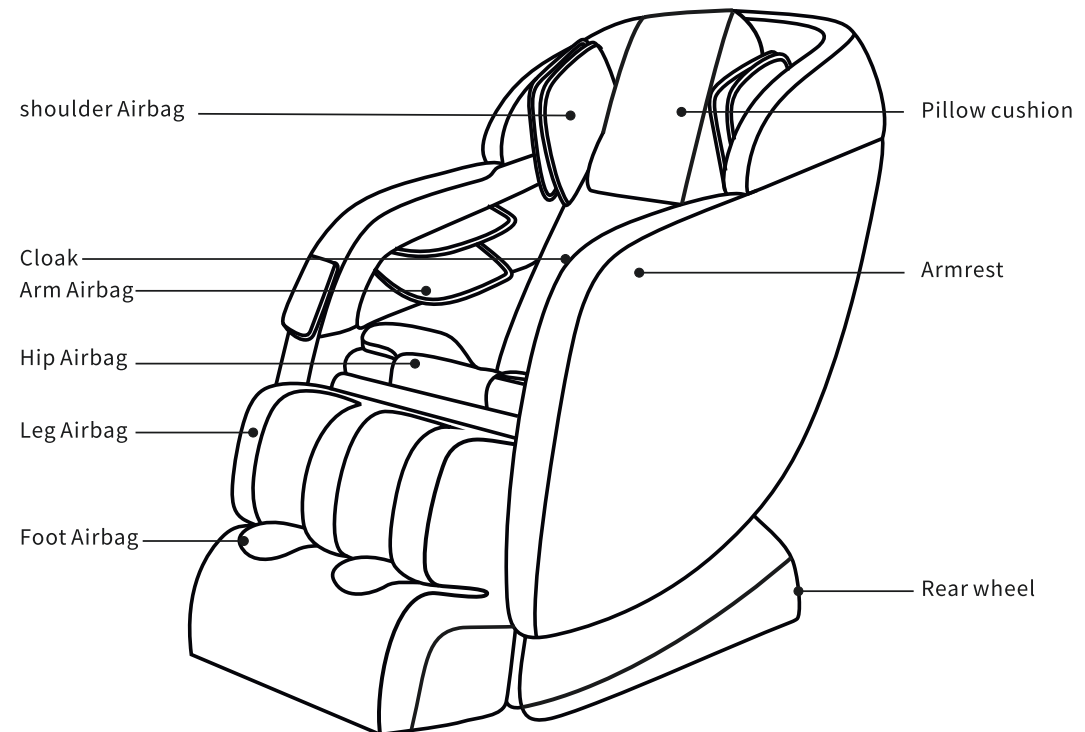
02 UNPACKING AND INSTALLATION

Remove the outer packaging carefully

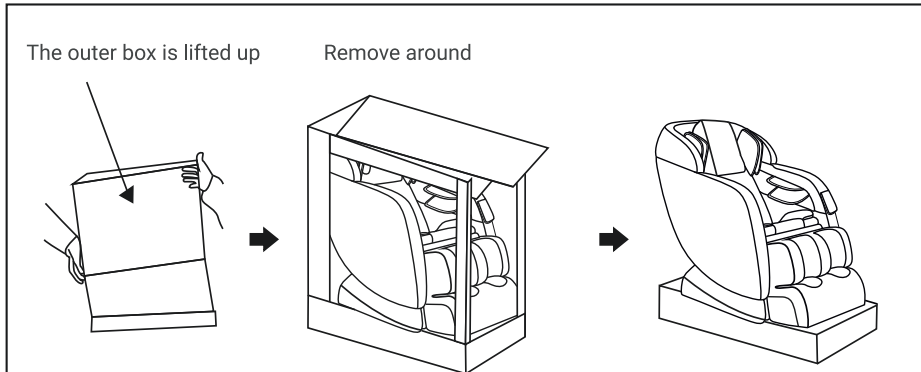Position the Massage Chair

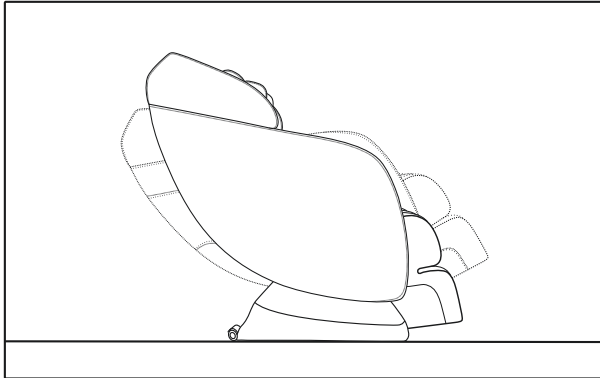- Carefully lift the massage chair out of the base carton and place in desired location
- When positioning the massage chair, be sure to leave 15cm space at the back of the chair and 15cm space in front of the chair so that the legs and backrest can be reclined and raised without contacting walls or other items
- Do not place the power cord in the electric control box at the rear of the chair. Check to ensure the cord is not under the massage chair.

03 INSTALLATION


Attach the power cord to the socket at the back of the Chair and plug into the mains.

ATTENTION:

- Do not use it in wet areas such as bathrooms or saunas.
- Place the chair on a flat, non-slippery surface.
- Do not place in direct sunlight or other hot or dusty locations


04 MOVING THE CHAIR

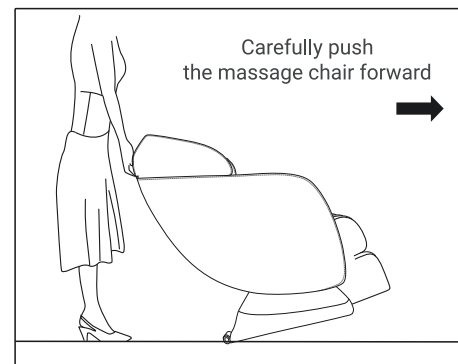- Place both hands on the back of the chair and apply weight to tilt to a 45-degree angle. Move the chair carefully on the roller castors to locate in the desired position.

- Note: Do NOT attempt to move the chair while someone is sitting in it.

05 REMOTE CONTROL MANUAL

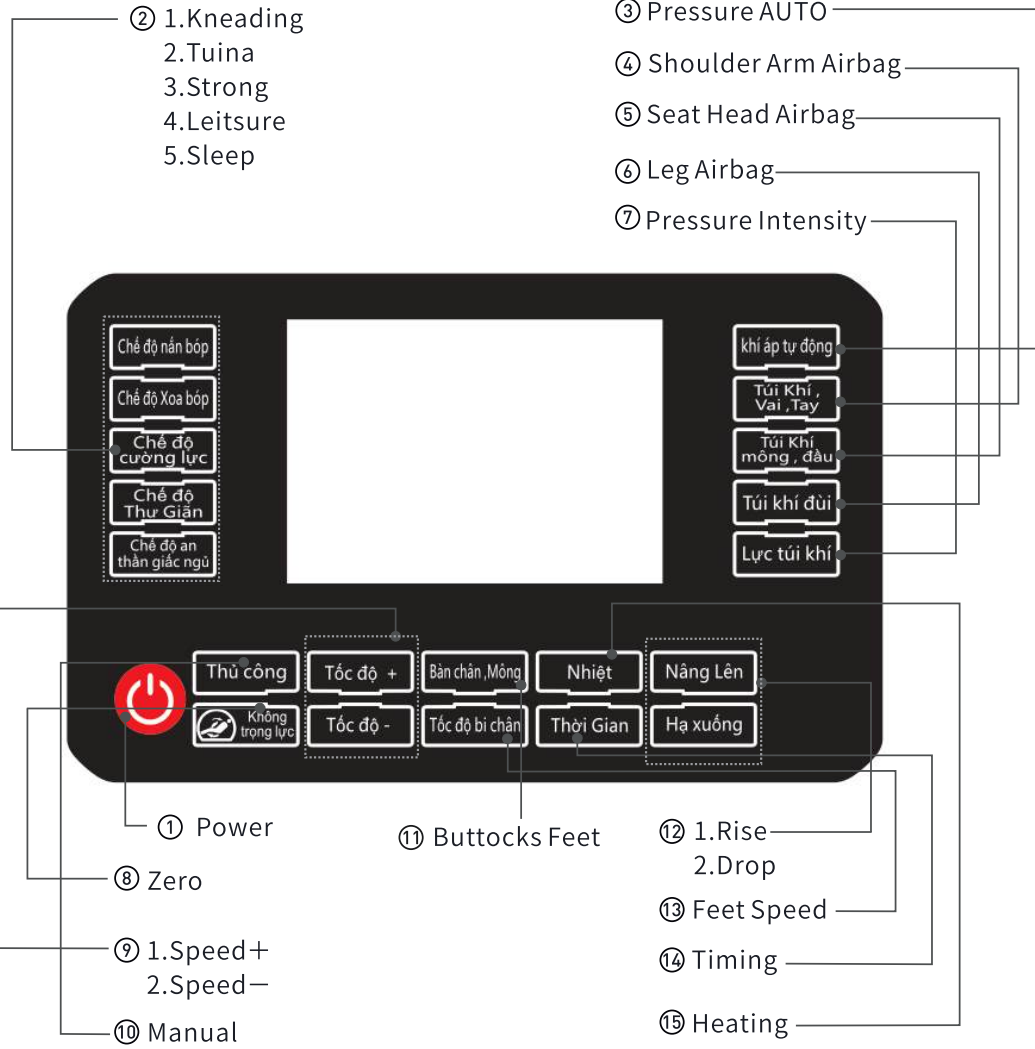06 REMOTE CONTROL KEYS

 **Power:** Press to turn the Massage Chair On or Off.

Press Power: Press the power button and select either automatic message or manual message option. The base time is set to 15 minutes.

Automatic Mode Program Selection:

-  **Kneading mode**
-  **Tuina mode**
-  **Strong mode**
-  **Leisure mode**
-  **Sleep mode**

-  **Pneumatic automatic:** Press to turn the pneumatic automatic On/Off.
-  **Shoulder arm balloon:** Press to turn the shoulder arm balloon On/Off.
-  **Seal cushion head air bag:** Press to turn the seal cushion head air bag On/Off.
-  **The leg bag:** Press to turn the leg bag On/Off.

 **Zero:** Press to move the chair into the zero-gravity position (3 steps)

 **Soles of the feet hip:** Press to turn the foot roller On/Off or adjust the speed (3 levels)

 **Foot speed:** Press to turn the foot speed on.

 **Timing:** Use this button to change the massage time anywhere between 5 and 30 minutes. The default time is 30 minutes.

 **Heating:** Press to turn the heating On/Off.

 **Rise:** Long press, calf rise; backrest down.

 **Drop:** Long press, lower leg; backrest rises.

08 PRODUCT SPECIFICATIONS

Product Name	Luxury Massage Chair
Model	MK-5150
Rated Voltage	AC220V – 60V
Rated Power	90W
Weight	65kg