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
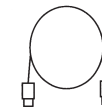
Cáp sạc

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Rảnh tay

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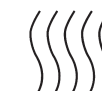INSTRUCTION MANUAL
NECK MASSAGER BELT
BUHEUNG **MK-332**

Scan the Zalo QR code to
activate online warranty

Thank you for choosing and supporting BUHEUNG KOREA!

We hope that you are satisfied with our product and service, and that we have the opportunity to serve and accompany you and your family.

If you have any feedback or reviews regarding the quality of our product and service, please contact our hotline at **1800 8379** so that we can receive and promptly address them.

In addition, please scan the **Zalo OA** code to activate the electronic warranty, receive additional gifts, and enjoy **promotions** and **24/7** support from BUHEUNG.

Thank you once again and we are happy to have take care of you and your family's health!

A stylized, cursive signature of the word "Buheung" in black ink.

TABLE OF CONTENT

- | | |
|-------------------------|-------------------------|
| 1. TROUBLE SHOOTING | 6. CHARGING METHOD |
| 2. ATTENTION | 7. PROHIBIT |
| 3. REMOVE THE BATTERY | 8. ATTENTION TO USE |
| 4. PRODUCT INFORMATION | 9. PRODUCT PROTECTION |
| 5. THE RIGHT WAY TO USE | 10. ESSENTIAL KNOWLEDGE |

1. TROUBLE SHOOTING

If a fault occurs during operation, perform the following checks. If the fault or the product cannot be removed. The product is still not ideal operation, please cut off the power supply, should be repaired by professional and technical personnel, do not disassemble This product!

Basic fault	Possible causes	Basic exclusion method
Unable to boot	Low charge	Please charge in time
Weak massage	Low charge	Please charge in time
Non-hot compress	Hot compress function is not enabled	Please turn on the hot compress function

2. ATTENTION

Before disposing of this product, the battery shall be removed as required by regulations;

When removing the battery, the product must be powered off;

Batteries shall be disposed of safely in accordance with laws and regulations;

Non-professionals please do not use yourself!

3. REMOVE THE BATTERY

1. Cut off the external power supply and remove the product shell;
2. Unplug the battery.
3. Remove the battery.
4. Waste rechargeable batteries should be placed in a fixed recycling point;
5. Please comply with local environmental protection laws and disposal guidelines to protect the environment.

The name and content of hazardous substances in this product are marked as follows

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This form is prepared in accordance with the provisions of SJ/T11364

Circuit board components include components of printed circuit board machines, such as resistors, capacitors, integrated circuits, Connectors, etc. 0: indicates that the content of the harmful substance in all homogeneous materials of the part is GB/T26572 Below the limits specified in 572. X: indicates that the hazardous substance is present in at least one homogeneous material of the part In excess of the GB/T26572 limit requirements.

ENVIRONMENTAL PROTECTION LIFE DESCRIPTION

This electronic and electrical product shall be used under the conditions of use described in the manual. contained Harmful substances or elements will not leak out, will not cause pollution to the environment or human, The time limit for causing serious damage to property is 10 years.

STORAGE

When the storage period exceeds 3 months, the machine should be recharged and used after unpacking. Place not Can be close to the ground, walls and roof. The room should be well ventilated to avoid strong sunlight and corrosive gas.

TRANSPORT

Product packaging with simple shock-proof facilities, suitable for air, rail, road and ship transport should be avoided Rain, snow splash, inversion and collision.

POWER CHECK

Check the power supply. If the power supply voltage exceeds the specified range, do not start the power supply

UPKEEP

Measures shall be taken to prevent moisture, rodents and pollution. Cleaning should be carried out in the state of power off. UsePlace the instrument in a safe place and do not place heavy objects on the instrument. When something goes wrong Professional and technical personnel should be asked to repair, non-professional personnel can not disassemble the product.

DISPOSAL AND RECYCLING

This product should not be discarded with household waste. When you decide to dispose of this product, do so in accordance with local environmental laws and guidelines.

WARRANTY CARD

Customer: Contact number:.....

Address:.....

Where to buy: Date of purchase:years month day

Product information: MODEL number: S/N:

This card is for service after purchasing massage series products, please keep it properly, and please refer to it in detail Fill in this information card for prompt service. If you need service, please contact original purchase Shop number or our company contact, and show the original purchase store sales list, when we are happy to serve you.

CERTIFICATE OF CONFORMITY


Product name: **Shoulder & neck massager** Model number:

Inspection number: Production date:

Declaration: With the technical innovation and improvement of the product, the same type of instructions as this product Subject to change and updated without notice. The legend shown in this manual is standard Based on the appearance of the model massager, it is for illustrative purposes only, the actual appearance is to be purchased The product shall prevail. Due to product improvement, some parameters may be changed, the changed parameters of this product Products have no retroactivity. The final interpretation of the relevant contents of this product manual belongs to the manufacturer.

4. PRODUCT INFORMATION

Product name	NECK MASSAGER BELT		
Product model	MK-332	Rated voltage	5V
Rated power	26W	Battery capacity	2000mAh

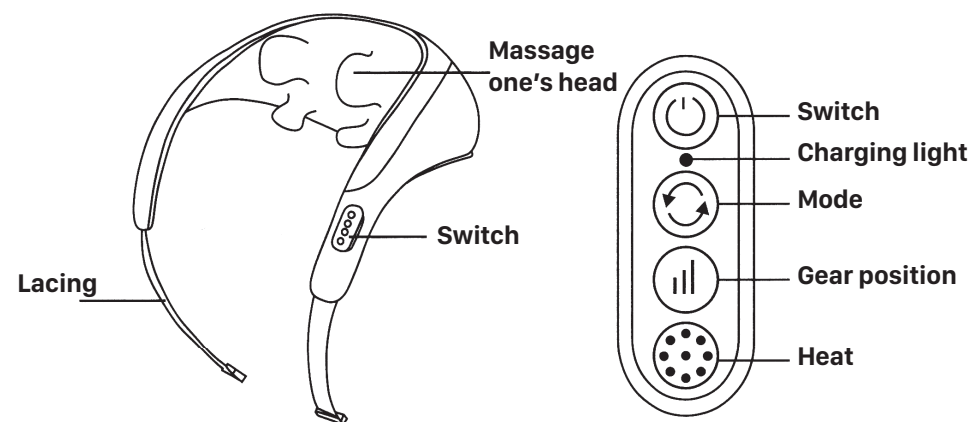5. THE RIGHT WAY TO USE

PREPARATION BEFORE USE

Please check whether the product and accessories are complete and damaged.
Press and hold the on key for 2 seconds, and the product starts to work.

* Mode selection

The default mode is 1. Press the mode key to switch between mode 1, Mode 2 and Mode 3.

* Strength mode

Boot the default second gear, short press the gear button, in turn for a one-one- three gear one-two gear cycle.


* Heating mode

Boot the default heating off, short press the heating key, heating high and heating low in turn and turn off the heating, so as to cycle.

* Product shutdown

Automatic shutdown: After 10 minutes of use, the product stops working

Manual shutdown: press the on key, the product stops working.

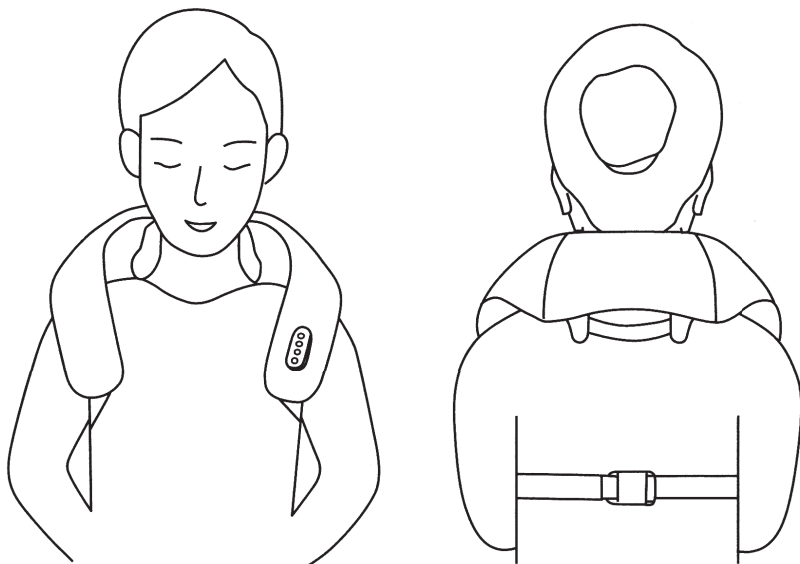
6. CHARGING METHOD

1. This product cannot be used while charging
2. Plug the charging cable into the charging port on the product.
3. When the product is working, it is connected to the power supply to charge, and the product starts the self-protection program Dynamic stop working.
4. The product charging indicator light is red, and after charging is completed, the indicator light is green.

* It is normal for the main and power adapter to heat up slightly during charging.

END OF CHARGE

After charging, unplug the power supply.

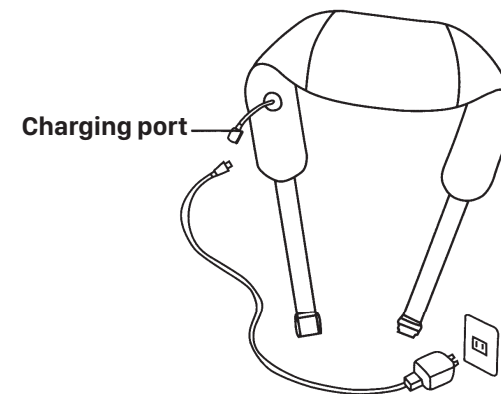7. PROHIBIT

1. Please read the instructions carefully and do not use the product continuously for a long time to prevent excessive pressing Rubbing causes muscle pain.
2. Do not use when lying flat or half lying.
3. Do not use this product during precise operation or driving.
4. The massage site has been operated on or has skin injury, redness, inflammation, congestion and other discomfort symptoms Do not use.
5. Bleeding in patients with heart disease or history, acute diseases, skin diseases, and malignant tumors People who are prone to blood diseases, osteoporosis, soft tissue diseases, stroke and other physical abnormalities please Do not use.
6. Do not use for people with slow movement, physical and mental dysfunction or intellectual disabilities. If necessary, please follow the instructions Follow the doctor's orders.

7. Pregnant women and menstrual period are not recommended to use massagers, if necessary, please follow the doctor's advice.

8. When storing the product, please keep it in a place that is not accessible to children. 9.8 Children under the age of 8 to 12 years old must be used in the hands of a parent or a person capable of supervision Use under guidance and care.

10. Do not scratch or stab this product with sharp objects.

8. ATTENTION TO USE

1. This product is a household electronic health care product, not as a medical device is not suitable for diagnosis, treatment and treatment Other commercial uses.

2. Before use, please remove accessories or other items that may affect the massage

3. Please stop using the massage if there is pain or discomfort.

4. When using the product, if there is an abnormal phenomenon, please immediately shut down, cut off the power and contact customer service Hotline, please do not disassemble and repair yourself.

5. There is a hot surface in the massage area, and people who are not sensitive to heat need to pay attention when using.

6. When used in other specific environments, please comply with safety regulations.

7. The appliance has a hot surface, which must be paid attention to when used by people who are not sensitive to heat.

8. This device is not intended to be caused by physical infirmity, sensory or mental impairment or lack of relevant experience The knowledge is used by people (including children) unless there is someone responsible for their safety who uses it for them This product provides supervision and guidance.

9. Children should be taken care of to ensure that they do not play with this product

9. PRODUCT PROTECTION

1. Please use the charging cable configured by our company and the adapter that meets the standard to charge the product. And pay attention to the safety of electricity.

2. Keep this product away from water, fire and corrosive environment.

10. ESSENTIAL KNOWLEDGE

* ABOUT FREQUENCY OF USE

1. Use as needed, within 10 minutes at a time. It is recommended to use once and take a rest before using.

2. If used excessively, it may cause physical discomfort and skin injury.

3. Therefore, when using for the first time, it is necessary to observe the physical condition and changes in the skin, and try to make it short With time.

*NOTE BEFORE USE

1. When using wood products, please adopt a personal comfortable position and completely relax your body.

2. If you feel pain or discomfort on the skin, stop using and seek medical attention immediately.

* ABOUT THE STATE AFTER USE

After a long massage, there may be redness and fever in the massage area while there is comfort This is a normal phenomenon.

* PRODUCT CLEANING

To clean the appearance, wipe the surface with a cloth dipped in water or 75% alcohol. Do not rinse with water.


PRODUCT CONFIGURATION

Host machine

Specification

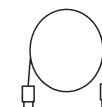Charging cable

PRODUCT CHARACTERISTICS

Bionic massage

Free hands

Light compress