

1234

A B C D E F



此页仅为材质尺寸说明
无需打印

- 技术要求:
- 1、印刷表面不能模糊不清、不得有漏印、错印、混色、重影、色差、色薄、花痕、表面异物等影响外观的不良现象。
 - 2、文字、图案印刷不得有断字、断线、残缺、粗细不均、字体粗细偏差、锯齿等不良现象。
 - 3、印刷表面不得有斑点、针孔、边角漏印、桔皮等不良现象。
 - 4、印刷需符合设计文档。
 - 5、环保要求:需符合RoHS指令
 - 6、具体内容参考电子文档。

																艾力斯特健康科技有限公司			

A B C D E F

CONTACT US

If your massage chair encounters any issues during use,
please contact us for support and warranty service.



🌐 WEBSITE: poongsankorea.vn

☎ HOTLINE: 1900 2296

MORE
25 *Year*
COMMITMENT
WITH VIETNAMESE FAMILIES



POONGSAN

MCP-279 MASSAGE CHAIR



INSTRUCTION MANUAL

Please read the instructions carefully before use.

POONGSAN

POONGSAN

Care without a day off



25^{Years}
COMMEMORATE

WITH VIETNAMESE FAMILIES

WARRANTY UP TO
5 YEAR

EASY RETURN
in
7 DAYS

0% INSTALLMENT
INTEREST RATE

GIFTS OF GRATITUDE RESERVED FOR
**LOYAL
CUSTOMERS**



poongsankorea.vn

HOTLINE: 1900 2296

By SCX

A JOURNEY FULL OF PRIDE AND EFFORT

Impressive Figures

With more than 20 years of development, Poongsan has been alongside consumers, achieving impressive milestones and affirming its efforts to provide an exceptional experience with high-quality massage chairs.

20+

20 Years of Partnership

Poongsan is a reputable brand from South Korea, first introduced to Vietnam in 2001.

68+

68 Showrooms Nationwide

With a network of showrooms across the country, Poongsan proudly stands as a trusted destination for direct shopping, offering convenient experiences and dedicated service to all customers.

99+

99% Customer Satisfaction

Poongsan is committed to outstanding quality and dedicated service, ensuring that every customer is satisfied with their experience.

15m

15 Million Customers

Poongsan is proud to be a leading and trusted brand, chosen by millions of customers, affirming its position as a top brand in the journey of health and family care across the country.

DEDICATED SERVICE

BRINGING THE BEST VALUE TO YOU



FREE SHIPPING
IN 24 HOURS



FREE
EXPERIENCE



LIFETIME
MAINTENANCE



DEDICATED
WARRANTY SERVICE



POONGSAN

We are committed to providing customers with comprehensive warranty services, including:

Product Warranty: All Poongsan massage chairs come with an extended warranty, ensuring the stable and effective operation of every component.

Technical Support: Poongsan's professional technicians are ready to assist customers in quickly and effectively resolving any technical issues.

Maintenance and Repair: Poongsan offers regular maintenance services and repairs as needed to ensure that the products continue to perform at their best.

NORTH REGION

Ha Noi	Manufacturer's Warranty Center	2020 Nguyen Son, Long Bien, Ha Noi
Hue Binh	Hue Binh authorized warranty Station - Tuan Cuong	151 An Quang Vuong, To 3
Lang Son	Lang Son authorized warranty Station - Hai Binh	369 Ba Thau Street
Lai Chau	Lai Chau authorized warranty Station	To 6, Ward Quyet Tien, Lai Chau City, Lai Chau Province
Son La	Son La authorized warranty Station - Thai Tung	To 2, To Hieu Ward
Dien Bien	Dien Bien authorized warranty Station - Tran Ngoc	88, To 6, Thanh Nam Ward, Dien Bien City
Tuyen Quang	Tuyen Quang authorized warranty Station - Tam Viet	174 Binh Thuan Street, To 10, Tan Quang Ward, Tuyen Quang City
Yen Bai	Manh Su Electronics - Refrigeration Co., Ltd.	52 Tien Hung Dao, Hong Phong Residential Area, Hong Ha Ward, Yen Bai City, Yen Bai Province
Bac Ninh	Bac Giang authorized warranty Station	243 Giep Hai, Dinh Ke Ward, Bac Giang City
Quang Ninh	Minh Thang authorized warranty Station	To 2, Khu 5, Ha Tu Ward, Ha Long City, Quang Ninh Province
Thai Nguyen	Quang Hoa authorized warranty Station	298 De Gang Thiep, To 2, Cam Gia Ward, Thai Nguyen City
Bac Giang	Bac Giang authorized warranty Station	243 Giep Hai, Dinh Ke Ward, Bac Giang City
Vinh Phuc	Truong Phat Electronics Co., Ltd.	103 Lam Son Street, Tich Son Ward, Vinh Yen City, Vinh Phuc Province
Hai Phong 1	Binh Minh Warranty Center	258 Mai Chay, Ngo Quyen District, Hai Phong City
Hai Duong	Trung Dat HO Trading & Service Co., Ltd.	8, Alley 78, Cong Cau Street, Hai Tien Ward, Hai Duong City, Hai Duong Province
Hung Yen	Quat Tien authorized warranty Station	233 Dien Binh Street, La Loi Ward, Hung Yen City
Nam Dinh	Thanh Tong Technical Service Co., Ltd.	825 Dien Bien Street, Loc Hoa Ward, Nam Dinh City
Ha Nam	Thanh Thang Warranty Station	QL 1A Le Huat Street, Minh Khai Ward, Phu Ly City, Ha Nam Province
Phu Tho	Tang Tu Refrigeration Co., Ltd.	1827 Hung Vuong Avenue, Gia Can Ward, Viet Tri City, Phu Tho Province
Ninh Binh	Phung Thinh	409 Hai Thuong Lam Ong, Le Loi Street, Van Giang Ward
Thai Binh 1	Binh Long Electronics & Refrigeration Co., Ltd.	15N Nguyen Tang Quat, Tran Lam Ward, Thai Binh
Thai Binh 2	Thai Binh Electronics Trading Co., Ltd.	La 27, Tran Thu Do Street, To 9, Tien Phong Ward, Thai Binh City
Hai Duong	Nguyen Dinh Huan	Phuong Di Village, Vinh Hung Commune, Binh Giang District, Hai Duong
Thanh Hoa	Minh Anh Trading & Service Co., Ltd.	543 Ba Thieu Street, Dong Tho Ward, Thanh Hoa City
Nghia An	Quat Vinh Refrigeration Business	7042 Hermann, Yen Vinh Block, Hung Phuc Ward, Vinh City, Nghia An
Ha Tinh	Huu Thinh Ha Tinh Trading & Service Co., Ltd.	10, Alley 581, Mai Thuc Loan Street, Dong Man Commune, Ha Tinh City

A Message of Gratitude

Dear Valued Customers,

First and foremost, we would like to express our sincere and heartfelt thanks to you for trusting, choosing, and using our massage chairs and other massage products. Your satisfaction and trust are invaluable sources of motivation, driving us to continually improve and provide high-quality products that support health care and enhance the quality of life.

We are proud to accompany you on your journey to find effective relaxation and health recovery solutions. Each of our massage chairs and massage tools is carefully designed with advanced technology, aimed at providing a wonderful relaxation experience and supporting health care in the best way possible.

In response to your trust, we commit to continually improving product quality, enhancing our services, and offering even greater value. Your satisfaction is our greatest goal, and we are always striving to achieve it.

Once again, we sincerely thank you for choosing our brand. We look forward to continuing to receive your trust and support in the future.

Warm regards,"

About the Brand

POONGSAN is a renowned brand from South Korea, present in Vietnam since 2001.

After more than 20 years, Poongsan has established itself as the top choice for many Vietnamese families, offering comprehensive health care solutions. Poongsan specializes in manufacturing and distributing medical devices, massage equipment, physical therapy products, and health care items.

Founded in 1985 as a South Korean corporation, Poongsan boasts over 35 years of experience and is proud to lead the way in applying modern technology to the health care industry. Each Poongsan massage chair is designed with advanced technology, recreating a professional massage experience akin to the hands of a skilled therapist.

Thanks for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay more attention to some important safety warning and keep this manual book properly for future use.

Note: Our company reserves the right to change and explain the design of the product. If there is some difference between the product function in market circulation and the description in the manual, it is a normal update of the product, and our company will not inform you otherwise. The specific product function and color shall prevail in kind.

TABLE OF CONTENTS

Before use

Safety Precautions	02-05
Installation parts	06
Place massage chair	07-08
Before use	09-10
Product structure	11-14

Usage method

Operation Instructions	15-20
Massage close	21

Maintain

Cleaning and Maintenance	22
--------------------------	----

Important

Troubleshooting	23
Product Specifications	24

Safety Precautions (Must be observed)

- Please read all relevant instructions carefully before use.
- The notes marked here will help you to use the product safely and correctly.
- In order to clearly indicate the extent of harm and damage, precautions are divided into **"Warning"** and **"Caution"** that may be caused by incorrect use. Both of them are important for safety and must be followed.

 Warning	Indicates that serious injury or death may result from improper use.
 Caution	Indicates that misuse may cause minor injury or property damage.

- Please pay attention to the symbol before each safety item.

 Forbidden	Absolutely Prohibited
--	-----------------------

The following populations should not use this product

 Forbidden	■ This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. ■ Children should be supervised to ensure that they do not play with the appliance. ■ The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance. ■ Patients with osteoporosis. ■ Patients with severe heart disease and wearing implantable medical devices such as pacemakers and defibrillators that are susceptible to electromagnetic waves. ■ People who have a fever (above 38 degrees). ■ Pregnant women or shortly after childbirth and during their menstruation period. ■ People who follow medical advice and are prohibited from receiving massage due to thrombosis, aneurysm, severe varicose veins, or other circulatory diseases. Within two weeks after the appearance of blood clots, it is strictly prohibited to massage any form of body parts that have blood clots. ■ People with injuries or skin conditions. ■ Unsupervised children under the age of 14, or anyone unconscious ■ People advised by their doctors to rest or those feeling unwell. ■ People with spinal abnormalities or scoliosis are prohibited from receiving this product due to medical advice due to other diseases. ■ People with wet bodies.
--	---

Safety Precautions (Must be observed)

Important safety warnings

 Warnings	■ If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard. ■ Keep children away from the movable parts of this product. ■ Use a properly grounded power source that is compatible with this device. ■ Unplug the power cord after use or before cleaning to prevent personal injury or damage to the device. ■ Use this device according to the instructions in the manual. ■ Do not use any parts or accessories that are not recommended. ■ Do not use this product outdoors. ■ Do not use near open flames such as stoves, and do not smoke while using. ■ Do not use heating equipment such as electric blankets on this product. ■ Do not sit on the back of the chair, leg massagers, shoulders, armrests, or stand on the seat. ■ Please read the instruction manual before operation. ■ Never use this product for any purpose not listed in the manual. ■ Limit each use to 20 minutes. ■ Do not use this product if the cover or leather is broken or damaged. ■ Do not drop anything into this product. ■ Never fall asleep while using this product. ■ Do not use this product if you are intoxicated or feeling unwell. ■ Do not use this product within one hour of eating. ■ Do not apply excessive force when using this product to avoid injury. ■ Do not use for purposes other than massage. ■ Precautions: Keep away from children. Do not store near high temperatures or open flames. ■ Fully insert the power plug into the outlet and ensure secure contact during use. ■ When using an extension cord, only use nationally certified products. ■ Disconnect power when leaving the premises. ■ Immediately discontinue use and disconnect the power if abnormal noise, odor, overheating, or malfunction occurs. ■ Do not power the unit from lamp socket outlets. ■ Total current draw of connected appliances must not exceed the rated current of the main outlet/circuit. ■ Never use near flammable/explosive materials. ■ Do not disassemble this product or modify its internal wiring, components, or safety features. ■ When unplugging, do not pull on the power cord. ■ Only use specified fuse ratings - never substitute with conductive materials or incorrect ratings. ■ Do not alter the trip settings of circuit breakers or residual current devices.
---	---

Safety Precautions (Must be observed)

Safety precautions	
 Caution	■ Check whether voltage is consistent with the specification of this product. ■ Never insert or unplug the power cord with wet hands. ■ Do not allow water to enter this product to avoid electric shock or damage. ■ Do not pull the power cord when plugging or unplugging. Avoid rough handling. ■ Do not damage the wire or modify the circuit of this product. ■ Do not clean electrical parts, such as the power switch or plug, with a wet cloth. ■ Stay away from this product during a power failure to avoid injury if power is suddenly restored. ■ Stop using this product immediately if it malfunctions, and contact your local supplier. ■ Stop using this product if you feel unwell, and consult your health care practitioner.

Product care and maintenance	
 Caution	■ This product should be maintained by the service center designated by the manufacturer only. No user should disassemble or maintain it without permission. ■ Do not forget to turn off the main power switch after usage. ■ Do not use this product if the power socket is loose. ■ If the product will not be used for an extended period, roll up the power cord and store it in a dry, dust-free environment. ■ Do not store this product at a high temperature or near open flame. Avoid prolonged exposure to direct sunshine. ■ Please clean this product with dry cloth. Never use thinner, benzene or alcohol. ■ The mechanism of this appliance is specially designed and does not require special maintenance. ■ Do not puncture this product with sharp objects. ■ Do not roll or pull this product on any uneven ground. Instead, it should be lifted before moving. ■ Please use this product intermittently. Do not keep it working continuously for a long period of time. ■ Allow intermittent operation - continuous long-term use is prohibited.

Safety Precautions (Must be observed)

Storage conditions	
 Caution	■ Storage room temperature: 5 °C to 40 °C. ■ Relative humidity (RH) between 20% and 80%; free of corrosive gases, and a well-ventilated room.


Operating environment	
 Caution	■ Do not use this product in a highly humid environment, such as a bathroom. ■ Stop using this product immediately when the ambient temperature changes sharply ■ Do not use this product immediately after a sudden change in ambient temperature. ■ Do not use this product in environments with heavy dust or corrosive gases. ■ Do not use this product in confined spaces or areas with poor ventilation.

Solutions to ordinary malfunctions	
 Caution	■ If the motor sounds during operation of this product, this is the normal operation sound. ■ If the controller is not working normally, check if the power plug and socket are securely connected, and ensure the power switch is turned on. ■ If the rated working time is over, this product will automatically activate the power-off button; If this product keeps operating for a long period of time, the temperature protector will stop this product automatically. The product can be used again after half an hour. ■ Weight limit: 120kg. Heavier users should purchase our reinforced models.

Installation parts





Controller bracket installation instructions

After the screw thread on the support knob of the controller rest assembly aligns to the screw thread position of the support mounting seat on the armrest of the massage chair, just turn clockwise and screw tightly for fixation.

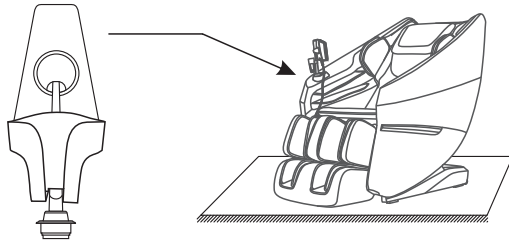

Calf installation instructions

1. After connecting the lower leg and the air tube on the massage chair with the wire harness, untie the fixed pipeline belt and insert the connecting part into the seat frame (Fig. 1).
2. Sleeve the lower leg fixing block on the round tubes on the left and right sides of the lower leg (Fig. 2).
3. Hang the lower leg sleeved with the fixing block on the hook of the massage chair seat frame (Fig. 3).
4. Fasten the lower leg on the massage chair with 4 M4 screws (2 on the left and 2 on the right). Check whether there is any abnormality in the up-and-down turning and whether (Fig. 4).

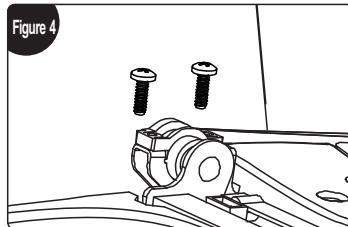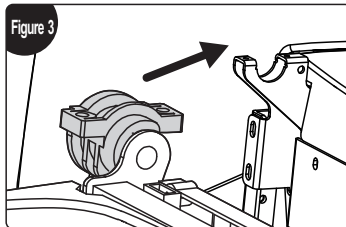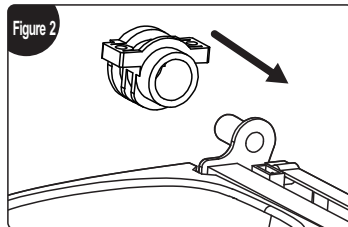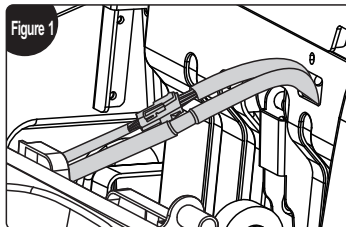


Place massage chair

Installation placement

Make sure there is enough space for the massage chair can be tilted at least 15 cm.



Do not expose the massage chair to direct sunlight or high-temperature environments (e.g., in front of heating devices). Place it in an area where it can avoid heat and direct sunlight. For long-term use in a fixed location, it is recommended to place a soft cushion underneath. The cushion should be large enough to cover the area where the chair and the footrest touch the floor.

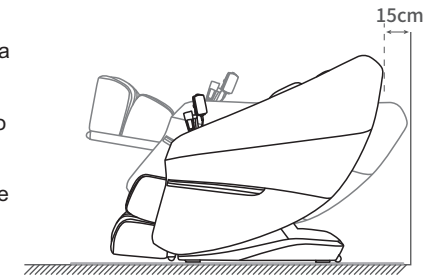


Product lift instruction

Two-person lifting required due to product weight.



1. Please lift the massage chair when moving it on delicate surfaces, such as wooden floors.
2. When moving on similar surfaces, lifting is recommended to avoid damaging the casters. Be cautious when lowering the chair into the designated position to avoid hand or foot injuries. Do not release the chair until it is completely flat. To prevent floor damage, it is recommended to place a mat or protective covering on the floor.

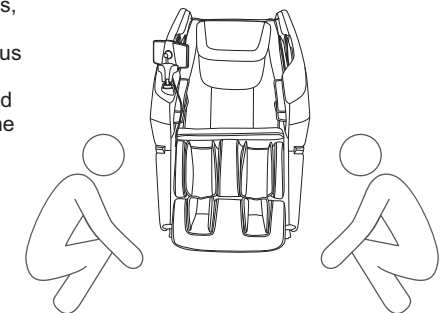


Place Massage Chair

Move method

Use the rear casters to move the massage chair. Tilt the chair backward approximately 45° as shown in the diagram to move it.



NOTE: Excessive force to tilt will cause the chair to tip over completely and product damage.

1. Before moving, ensure the massage chair is properly powered off, with both the backrest and calf sections fully reset, and the power supply disconnected.
2. The chair should remain empty during movement; no people, pets, or objects should be on it.
3. Do not move the chair over obstacles higher than 2mm or gaps wider than 5mm.
4. When moving the chair using the casters, avoid continuous movement over distances greater than 50 meters. If necessary, move the chair slowly and at a steady pace.

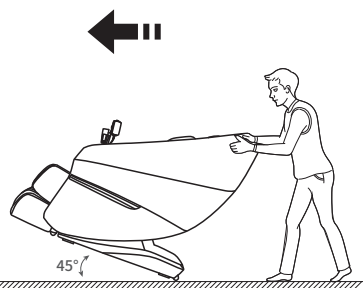

How to disassemble and use headrest cushion and backrest cushion


Using a headrest cushion can reduce the intensity of kneading massage on the neck and shoulders, and you can decide whether to use a headrest cushion according to your needs (recommended to use a headrest cushion). The backrest cushion and the backrest are connected by a zipper (1), and the headrest cushion and the backrest cushion are connected by a zipper (2).

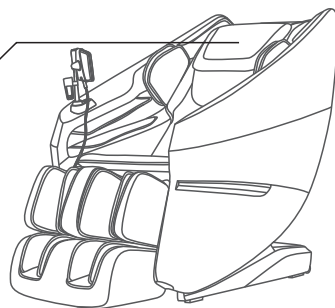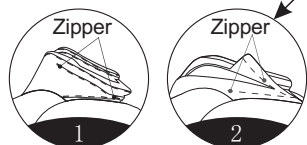

Before Use


Check around

Make sure that there are no people, pets and other objects around the unit.

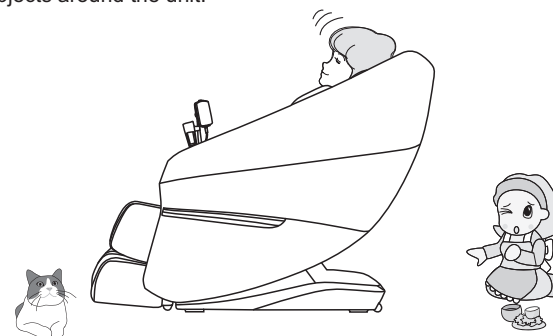

Make sure there are no foreign objects in the gaps of the machine

Before massaging, make sure that no foreign objects are caught in the legs, feet, etc. Under no circumstances allow hands, feet, and head to penetrate into the area between the lower leg mechanism and main body, or the area between the armrest and the capsule. Under no circumstances should hands, feet, or head sink under the fabric or leather covering of the massage area.

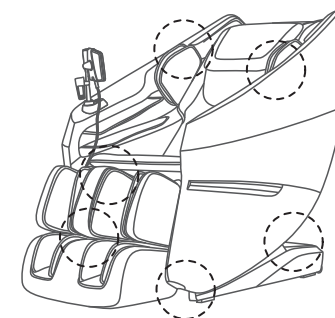


Before Use


Turn on the power

CAUTION


Plug in power cord

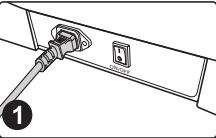

Schematic diagram of power on

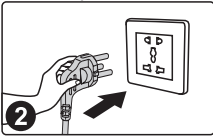

Turn on the power switch ("I" position)

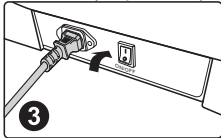


Warnings

- Please check the power cord and plug for damage before use.
- Please supply power to this product in accordance with the working power requirements of this product's calibration.
- It is strictly forbidden to use a power supply other than the nominal power supply of this product.
- It is strictly forbidden to connect various power conversion devices without authorization to supply power to this product.


Product Structure

Parts description

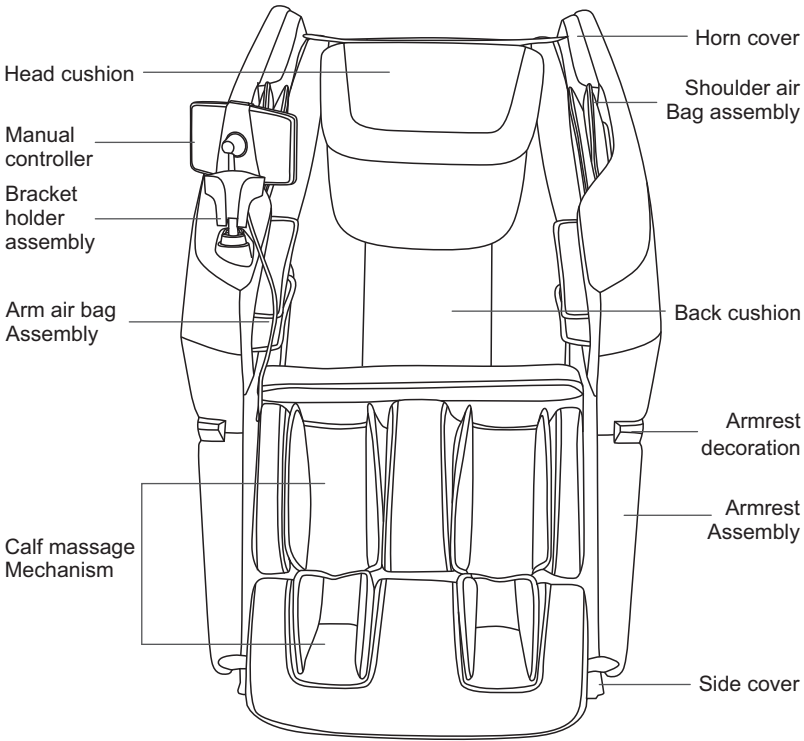

Product Structure


Function distribution

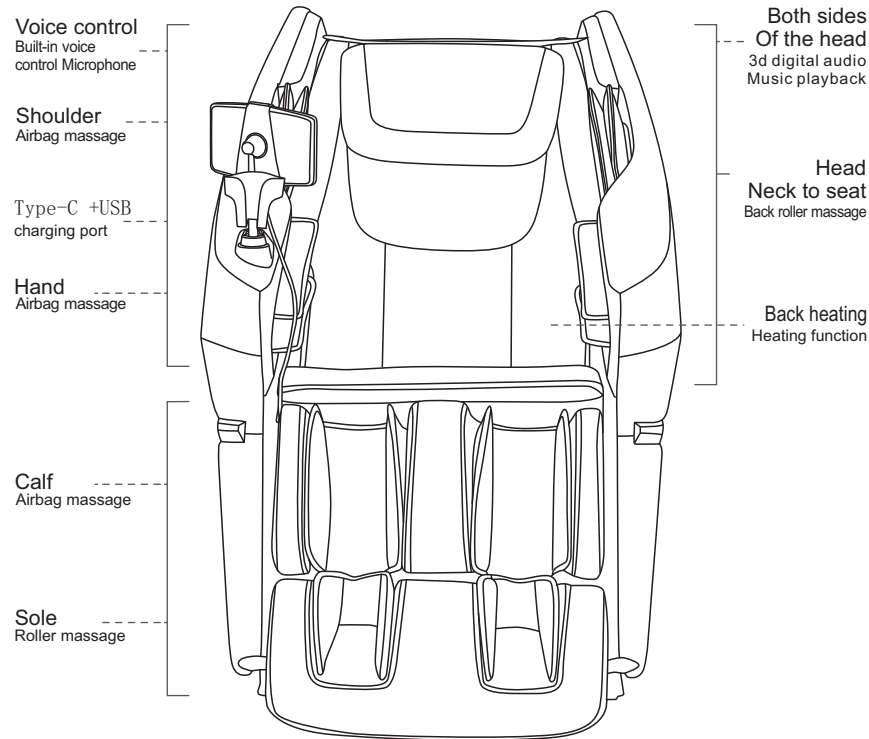

Product Structure


Introduction to the remote controller

1. When touching the manual controller screen keys in standby mode, the screen will light up.
2. After light up, you can only press the Lift/Lower Leg key and Get Up/Lying Down key to adjust the desired angle.

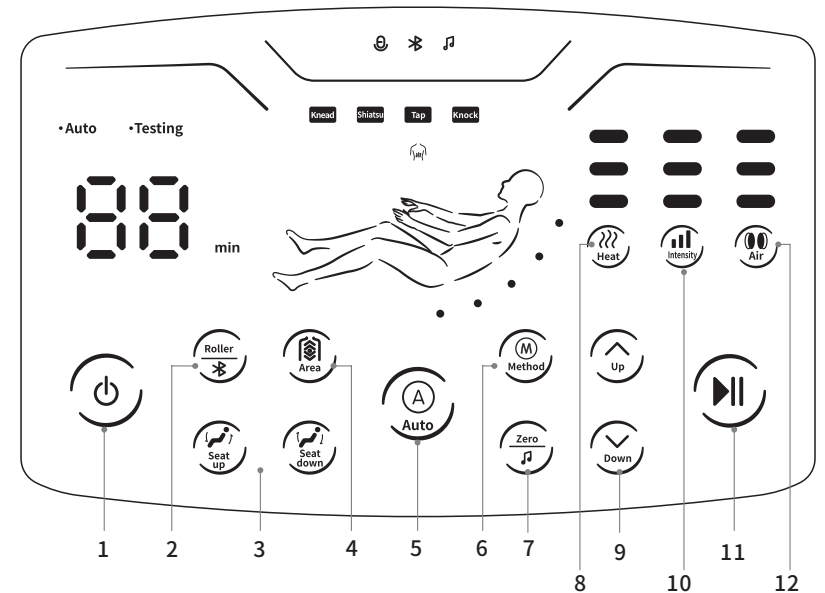

1. Standby key: Turn on or off the massage chair.
2. Foot roller: Short press turn on or off the foot roller massage function. Long press turn on or off the bluetooth.
3. Seat-up: Adjust the backrest frame up and the calf frame down.
Seat-down: Adjust the backrest frame down and the calf frame up.
4. Massage area: Adjust the massage area.
5. Automatic mode: Cyclic switching automatic mode.
6. Massage Method: Short press adjust the massage method. Long press adjust the volume.
7. ZERO: Short press: There are 3 zero-gravity modes and zero-gravity off, which can be cycled and switched as needed. Long press turn on or off the background music.
8. Heat: with a total of 3 levels that can be adjusted, or turn off the heating function.
9. Shoulder position adjustment: Manually adjust the position of human shoulder and massage the Shiatsu point.
10. Intensity: massage intensity adjustment, adjustable in 3 levels.
11. Pause: Short press pause or start the massage function. Long press return the massage chair to the sitting angle.
12. Air Intensity: Short press cycle adjustment of air pressure intensity, adjustable in 3 levels, turn off air pressure. Long press adjust the air part.

Product Structure

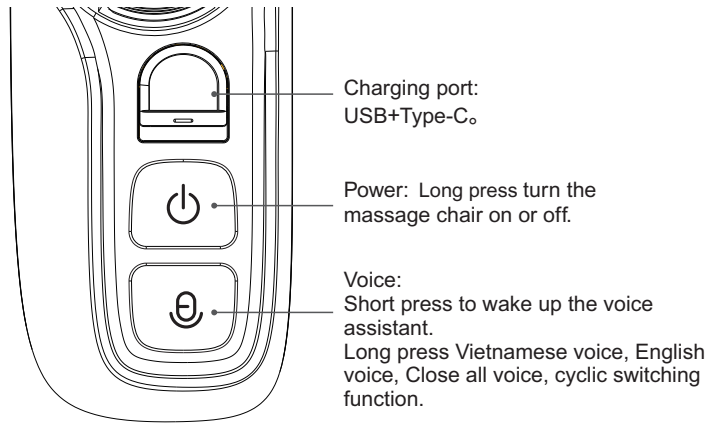


Operation Instructions

I. Start massage

- ①.Press the standby key  or long press the right armrest standby key  to start the massage.
- ②.Automatically adjust the angle.
- ③.Automatically detect human shoulder position, keep sitting posture fit with massage chair.
- ④.Conduct fine adjustment of shoulder position. If the detected position deviates from the actual shoulder position, adjust to the proper position by Up and Down keys for shoulder position. Shoulder position is adjustable. If no shoulder position adjustment within 10 seconds, the massage function will be started.

Note: During the detection process, you can choose automatic massage function keys to skip.

Button	Description
 Up	Adjust the massage ball position upward so that it reaches the proper shoulder position.
 Down	Adjust the massage ball position downward so that it reaches the proper shoulder position.

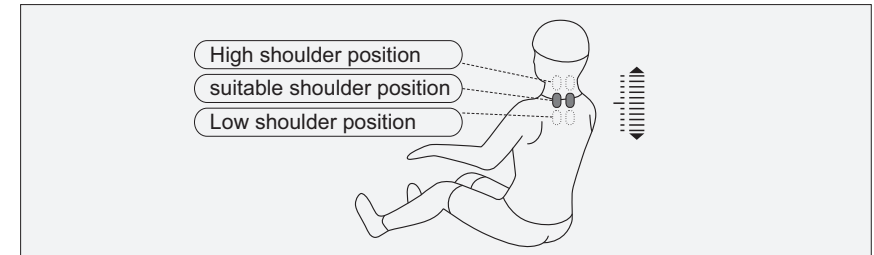

In the process of shoulder position detection, the "beep" sound of the manual controller is heard. At the same time, press the up and down keys of the shoulder position to adjust the shoulder position. At this time, you can also press automatic massage function keys to enter this mode for massage (Step 4 can be skipped directly).

• Pause

Function	Button	Description	Display
Pause		Short press when the massage function is running, press the "Pause" key to temporarily stop all massage functions; press the "Pause" key again to continue the massage function before the pause.	
		Long press return the massage chair to the sitting angle.	/

Operation Instructions

Auto Massage

Function	Button	Description	Display
Auto Massage	 Auto	Comfort Massage: The most comfortable massage experience. (With Background Music)	Auto 01
		Relax Massage: Gentle and comfortable deep massage, which can relieve body fatigue and relax the body. (With Background Music)	Auto 02
		Full body Massage: Deep massage, which can revive body vitality.	Auto 03
		Neck& shoulder Massage: According to the unique neck and shoulder massage method, through effective massage of cervical spine and scapula, neck and shoulder pain can be relieved and neck and shoulder fatigue can be improved.	Auto 04
		Stretch Massage: The lower leg frame and the backrest frame are lifted, effectively relieving tiredness, relieving joint soreness and restoring vitality.	Auto 05
		Back&Waist Massage: Focus on the back and waist to relieve back and waist soreness.	Auto 06
		Good-night Sleep: With kneading, shiatsu massage techniques, through the 12 health acupoints massage, help to regulate and balance on human body; so that people can fall asleep more quickly, improve sleep quality.(With Background Music)	Auto 07
		Sport Recovery: Quick recovery of body fatigue after exercising.	Auto 08
		Body Balance:Mainly massage the shoulders & neck, back & waist with kneading, kneading and shiatsu. Suitable for people who are mentally stressed, physically and mentally unbalanced. The massage intensity is gentle.	Auto 09
		Energy Boost: Deeply massage the whole body by means of kneading, tapping and Shiatsu techniques, to wake up the vitality of the human body and stimulate the energy.	Auto 10
		Fatigue Recovery: Restore the full body from fatigue and pain; achieve relaxation, improve microcirculation, and quickly relieve fatigue for the back and limbs.	Auto 11
		Hips&legs Fit: Massage effectively for buttocks and thighs toshape the body and reduce the accumulation of excess body fat.	Auto 12
		Spine Care: Protect the health of the whole spine, and improve the lack of exercise.	Auto 13

Operation Instructions

Auto Massage	 Auto	Lunch Break: Deeply massage the whole body by means of kneading, patting and finger pressing techniques, to promote blood circulation throughout the body and prevent blood clots. It is suitable for people who live at home for a long time and do not like to go out, and the massage intensity is moderately gentle.	Auto 14
		Sedentary: Aim to office sedentary group, increase the amount of physical exercise by massage, improve the body function.	Auto 15

Manual massage

Function	Button	Description	Display
Massage Method	 Method	Stop	/
		Knead	
		Tap	
		Knead and Tap	 
		Shiatsu	
		Knock	
Intensity	 Intensity	Adjustable 3 gears in intensity.	  
Massage Area	 Area	Fixed point: Massage in a fixed position.	H0
		Partial: Massage in a small area forth and back.	H1
		Upper back: Massage the upper back forth and back.	H2
		Lower back: Massage the lower back forth and back.	H3
		Full: massage the full body forth and back.	H4
	 Up	Massage position adjustment (in Fixed-Point and Partial mode, massage position can be adjusted).	Robot Up
	 Down		Robot Down

Operation Instructions

•Foot roller massage/

Function	Button	Description
Foot Roller	Roller 	Short press this key to turn on / off sole roller massage and it is adjustable of three gears in intensity, namely weak, medium and strong, it can be switched cyclically among three gears and off states. Long press turn on or off the bluetooth.

•Heat

Function	Button	Description	Display
Heat	 Heat	With a total of 3 levels that can be adjusted, or turn off the heating function.	

•Pneumatic massage

Function	Button	Description	Display
Airbag	 Air	Short press with a total of 3 levels that can be adjusted, or turn off the airbag pressure.	
		Long press adjust the air part. Full body airbag pressure (1), Upper body airbag pressure(2), Lower body airbag pressure(3), turn off the airbag pressure.	1,2, 3,OF

•Pneumatic massage

Function	Button	Description
Seat-up	 Seat up	Press and hold this key, the backrest frame rises slowly, the calf frame descends slowly, accompanying the sound of "beep... beep..." until reaching the top, if the key released the frames stop.
Seat-down	 Seat down	Press and hold this key, the backrest frame moves down slowly, the calf frame moves up slowly, accompanying the sound of "beep... beep..." until reaching the bottom. If the key released the frames stop.
Auto Seat-down	Zero 	Short press: There are 3 zero-gravity modes and zero-gravity off, which can be cycled and switched as needed. Long press: Turn on or off the background music.

•Other

Function	Button	Description
volume	 Method	Long press: adjusted the volume.

Operation Instructions

II. Voice Control Mode

- When the massage chair is turned on or in standby mode, the sound "Em ơi" will be produced to enable the voice control function.
- After such function is enabled successfully, the massage chair will reply "Dạ thưa chủ nhân." At this moment, once a corresponding voice control entry to be implemented is given within six seconds, the corresponding voice control function can be performed after a reply is received from the massage chair; if no correct voice control entry is given within six seconds, the voice control function automatically will be off automatically and needs to be enabled again.
- Note: When waking up the voice control function, please use English, the speech rate is moderate, Voice commands should be follow up the Voice entry.

Voice Control Entry		
English Voice Control Awakening Word: "Poongsan"		
No.	Voice Control Entry	Answer Entry
1	Massage on	Answer: Ok, massage on
2	Massage close	Answer: Ok, massage close
3	Comfort massage	Answer: Ok, comfort massage
4	Relax massage	Answer: Ok, relax massage
5	Full body massage	Answer: Ok, full body massage
6	Neck and shoulder massage	Answer: Ok, neck and shoulder massage
7	Stretch massage	Answer: Ok, stretch massage
8	Back and waist massage	Answer: Ok, back and waist massage
9	Open the air pressure	Answer: Ok, open the air pressure
10	Close the air pressure	Answer: Ok, close the air pressure
11	Up the seat position	Answer: Ok, up the seat position
12	Down the seat position	Answer: Ok, down the seat position
13	Change the other mode	Answer: Ok, change the other mode
14	Go little down	Answer: Ok, go little down
15	Go little up	Answer: Ok, go little up
16	Speed up	Answer: Ok, Speed up
17	Speed down	Answer: Ok, Speed down

Note: Command 14,15,16,17 can only be used after command 13.

Operation Instructions

Note: When waking up the voice control function, please use Vietnamese, the speech rate is moderate, Voice commands should be follow up the Voice entry.

Voice Control Entry		
Vietnamese Voice Control Awakening Word: "Em ơi"		
No.	Voice Control Entry	Answer Entry
1	Khởi động	Dạ Poongsan xin nghe ạ
2	Tắt máy	Vâng ạ, tắt máy
3	Nhẹ nhàng	Vâng ạ, nhẹ nhàng
4	Thư giãn	Vâng ạ, thư giãn
5	Toàn thân	Vâng ạ, toàn thân
6	Cổ vai	Vâng ạ, cổ vai
7	Kéo giãn	Vâng ạ, kéo giãn
8	Lưng eo	Vâng ạ, lưng eo
9	Túi khí	Vâng ạ, túi khí
10	Tắt túi khí	Vâng ạ, tắt túi khí
11	Nâng ghế lên	Vâng ạ, nâng lên
12	Hạ ghế xuống	Vâng ạ, hạ xuống
13	Đổi vị trí	Vâng ạ
14	Xuống dưới	Vâng ạ, xuống dưới
15	Lên trên	Vâng ạ, lên trên
16	Nhanh lên em	Vâng ạ
17	Chậm thôi em	Vâng ạ

Note: Command 14,15,16,17 can only be used after command 13

Operation Instructions


III. Operation Instructions for Bluetooth Speaker Play

1. Enter the smart device main menu. Open "Settings" and select the "Bluetooth" option. Open "Bluetooth", and Bluetooth will automatically "Search for a New Bluetooth Device".
2. The smart device will search for the corresponding Bluetooth model (e.g. MCP-279-XXXX). Click the "Match" button. If you are asked whether to match with the device, you just select "Match" to complete matching.
3. Turn on "Music Player" in the smart device and select the song to play.

IV. Massage close

During the massage, press the hand control standby key "⏻" or long press the right armrest standby key "⏻" to turn off. Will immediately turn off all massage functions, reset the backrest and leg, or when the massage time is up, the massage automatically turns off all massage functions, and the backrest mechanism and leg mechanism do not reset.

Cut off the power supply of the whole machine, as shown in the attached picture (the schematic diagram of the power cut off of the whole machine).

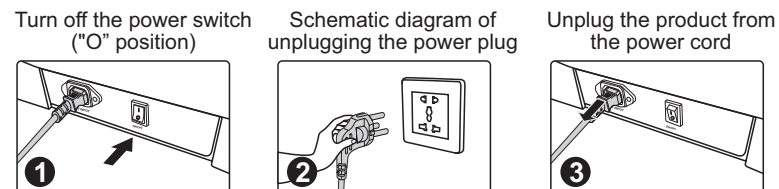

Cleaning And Maintenance

Cleaning synthetic leather and plastic parts

Wipe the leather with a soft, dry cloth.

- Do not use chemicals such as thinner, gasoline, alcohol, etc. to clean the unit.
- When using commercially available leather care products (rags), please follow the instructions for use.
- If synthetic leather is particularly dirty, wipe it in the following way.
 - ① Immerse a soft cloth in water or hot water containing 3% to 5% neutral detergent, then wring it out thoroughly.
 - ② Use the above soft cloth to wipe the leather surface.
 - ③ Wipe off the cleaning solution with a cloth soaked in clean water and thoroughly wrung out.
 - ④ Wipe with a wrung soft cloth.
 - ⑤ Let it air dry naturally.
- When it is difficult to remove stains, please immerse commercially available melamine foam in neutral detergent, and then wipe the machine with it.
- Do not use a hair dryer to dry the surface.
- Be careful when wiping the machine with denim and colored fabrics, as the color of the fabric will stain the surface of the synthetic leather.
- Do not keep these parts in contact with plastic for a long time, as this will cause discoloration.
- Synthetic leather may be discolored, so when using hair dye, please cover the touched part with a towel.



Troubleshooting

Fault	When the calf mechanism or the backrest reaches a certain position, press the corresponding adjustment button and it will no longer change, and a continuous "beep, beep" sound will be emitted.
Solution	Due to product structure design and ergonomic requirements, this product has a limited protection design for the limit state of attitude adjustment. At the same time, a reminder sound will sound when the limit is reached.

Fault	The function does not work or some functions do not work after the machine is turned on.
Solution	After the machine is turned on by pressing the switch key, it is necessary to manually select an automatic massage program or other working state. If it is not operated, the machine will automatically shut down after 20 minutes. If the manual operation mode is selected, the functions that are not selected may be disabled, and the working mode needs to be manually set in the manual mode.

Fault	The machine makes a leather rubbing sound when adjusting its posture.
Solution	This product adopts high-grade anti-corrosion leather that meets environmental protection requirements. Due to the inherent characteristics of leather materials, there will be a normal slight noise when rubbing at the relative movement, which is a normal phenomenon. Note: If you hear fabric or leather tearing sound when the machine is running, please stop using it immediately and seek professional maintenance.

Fault	The massage wheel did not reach the shoulders or neck.
Solution	If the head does not touch the head cushion or the back does not touch the backrest, the shoulder position may be lower than the actual position during the body shape detection process. Sit in the deepest position of the seat, rest your head on the headrest, and start over from the beginning.

Fault	The unit is damaged. The power cord or power plug is abnormally hot.
Solution	To prevent accidents, be sure to contact an authorized service center.

Please feel free to contact the designated service personnel for any questions or doubts about this product.

Product Specification

Model: MCP-279

Name: Massage Chair

Rated voltage: 220-240V~ 50Hz/60Hz

Rated power: 120W

Rated time: 20min

Safety structure: Class I

CONSULTATION & MANUAL USER



HOTLINE



1900.2296

*Please beware of scams
when making calls through unofficial channels.*