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Model: MK-321

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EYE MASSAGER

Model: MK-321

15 minutes a day, say goodbye to the eye fatigue!

**PRODUCTS ARE IMPORTED AND DISTRIBUTED EXCLUSIVELY BY TNHH
TM XNK NGUON SONG VIET STORE**

Ho Chi Minh : 235-237 Road No.9A Trung Son residential area,
Binh Hung commune, Binh Chanh district. Ho Chi Minh City.
Hotline : 18006977 | Service Centre : 028 7106 5858
Showroom Ha Noi (1st floor) : 2, Lane 382, Pham Van Dong , Co Nhue 2 ward,
Bac Tu Liem district, Ha Noi. Hotline : 024 37855966

- HOT COMPRESS
- ONE KEY CONTROL
- AIR PRESSURE
- BLUETOOTH
- 180 DEGREE FOLD
- LESS NOISE

1. BRIEF INTRODUCTION


Our products are with 180 degree folding design so it can fit all facial forms. Your eyes can be slightly pressed and massaged through the reasonable coordination between the air inflation and air deflation. In the meantime, your eyes shall be comprehensively and gently massaged through the hot compress

2. FIT FOR PEOPLE

1. With visual fatigue and pseudo-myopia.
2. Use eyes for long time with asthenopia.
3. Having been wearing glasses for long time.
4. Facing computer, cellphone or pad screen every day

3. NOTES FOR APPLICATION

1. The products are not fit for people with eye disease. The patients having been through the ocular operation or with ocular disease including glaucoma, cataracta and retinal detachment shall be prohibited for applying this product.
2. Patients with heavy stroke, hypertension or hypotension shall not apply this products before consulting the professional doctors
3. The products will be heated on the surface, people with less sensitivity of temperature must be aware of the heating surface of the products.

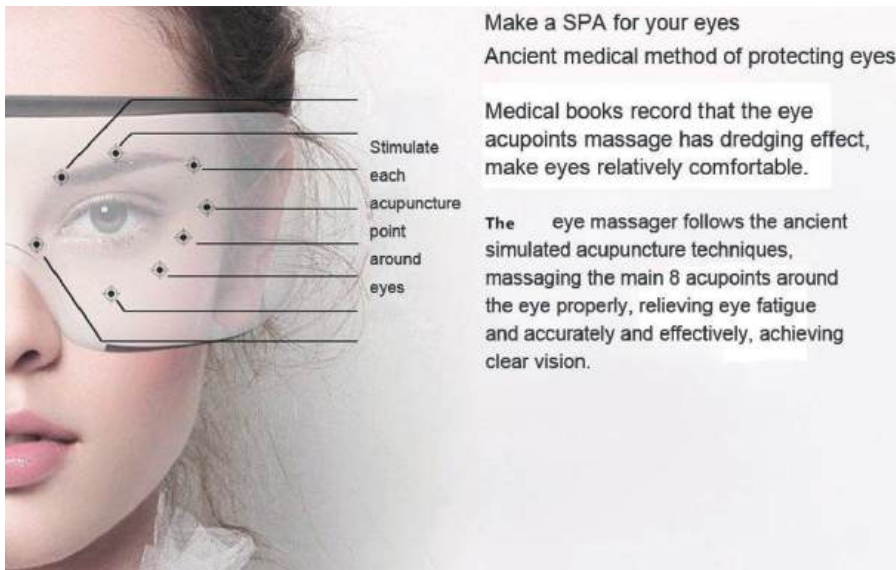OPERATING INSTRUCTIONS OF 4S

1. HOW TO USE?

Pressing the power button for 3 seconds and the product is switched on, the working status is hot compress; Pressing it again it will switch to air pressure status. Press the third time it will switch to air press and hot compress together. In this case air press and hot compress work at the same time. Press the power button the fourth time for 3-4 seconds it will be turned off. You can do this anytime you want.

2. HOW TO CONNECT IT TO BLUETOOTH?

Turned on the power and bluetooth signal icon will automatically blink, connect mobile phone bluetooth to **EYE MASSAGER**; Play music in mobile phone, adjust the appropriate volume, and enjoy it. Massagers can only connect one device at the same time, so if you want to switch device please turn off the device previously connected. Next time when you turn on it and the bluetooth will connect it automatically.

3. HOW TO CHARGE?

Charge through the android USB charging standard, the power indicator light flashing, after 2 hours fully charged, the light will stop flashing and normally on, please choose normal 5v charger, mobile phone charger or mobile power supply.

OPERATING INSTRUCTIONS

1. Pressing the power button for 3 seconds and the product is switched on. The working status is hot compress (Low temperature); Pressing it again it will switch to high temperature status.

2. Double click on the power button and it will switch to air press mode (low pressure). Then single click the button it will be changed to high pressure mode.

3. Double click on the power button again, it will be air press and hot compress (low) working at the same time. Single click on the button, high mode works.

4. Press the power button for 3-4 seconds it will be turned off. You can do this anytime you want.

5. How to connect it to bluetooth of your phone. Turn on the power and bluetooth signal icon will automatically blink, connect mobile phone bluetooth to SHU SEN; Play music in mobile phone, adjust the appropriate volume, and enjoy it. Massagers can only connect one device at the same time, so if you want to switch device please turn off the device previously connected. Next time when you turn on it and the bluetooth will connect it automatically.

6. How to charge

Charge through the android USB charging standard, the power indicator light flashing, after 2 hours fully charged, the light will stop flashing and normally on, please choose normal 5v charger, mobile phone charger or mobile power supply.

PROBLEM AND HANDLING METHODS

PROBLEM AND HANDLING METHODS	WHY	HOW TO DEAL
No indicating light	Low battery	Charge it for a few hours
No indicating light	Low battery	Charge it
Unable to power the products	Low battery	Charge it
Auto power-off after starting up	Low battery	Prompt changing
Unable to charge the product	Damaged usb cable	Change a cable
Unable to charge the product	Broken interface	Contact seller
Can not connect it to bluetooth of your phone	Turn off the bluetooth	Restart the massager, and connect it again

Learn more about Buheung at the website

www.buheung.vn

Read the manual carefully before use please!