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# INSTRUCTION MANUAL NECK MASSAGE PILLOW

BUHEUNG **MK-333**



Scan the Zalo QR code to  
activate online warranty

## PRODUCT SPECIFICATIONS

<b>Product Name</b>	NECK MASSAGE PILLOW	<b>Model</b>	MK-333
<b>Power</b>	10W	<b>Battery</b>	2000 mAh
<b>Operating Voltage</b>	3.7V (MAXIMUM 4.2V)	<b>Charging Voltage</b>	5V - 2A
<b>Charging time</b>	3.5 - 4 hours	<b>Usage time</b>	2 hours
<b>Product Dimensions</b>	245x230x105mm	<b>Weight</b>	0.8 kg

## PRODUCT STRUCTURE
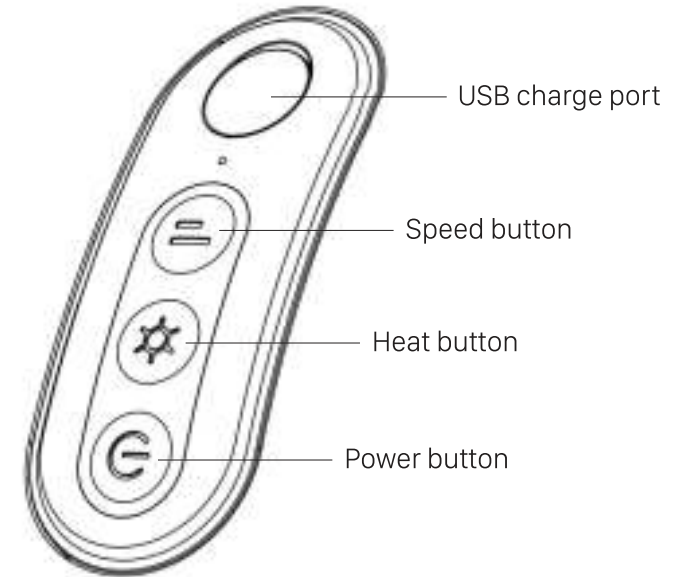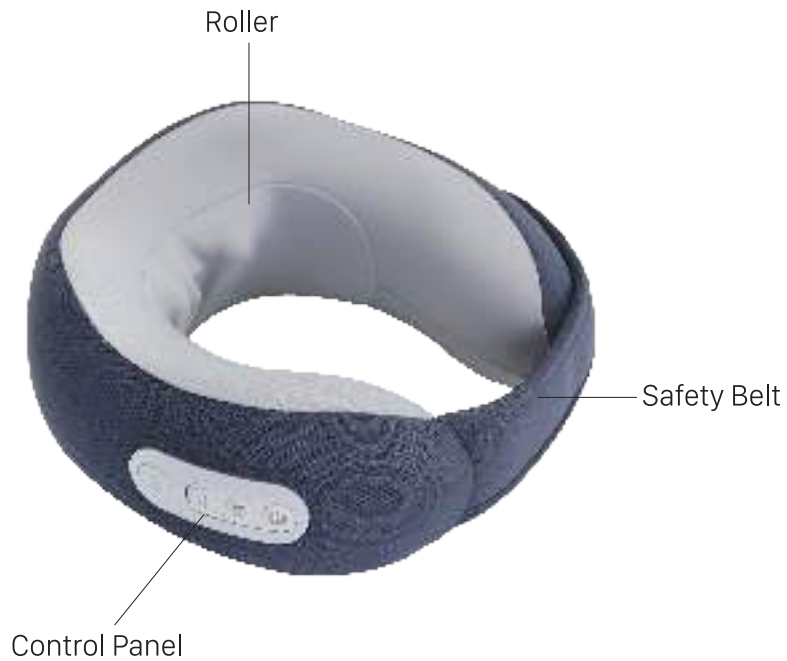
## INSTRUCTIONS FOR USE

Wear the RECHARGEABLE NECK MASSAGE PILLOW around your neck.

**Power On/Off:** Press the power button for 2 seconds to turn on or off the massage device.

Press the **Heating button** to turn on/off the heating function.

Press the **Speed button** to increase or decrease the massage speed.

The device will automatically change the massage direction every 1 minute.

The device will automatically turn off after 10 minutes of operation.

## CLEANING AND MAINTENANCE

1/ Do not use harsh chemicals or abrasive products to clean the device.

2/ Regular cleaning of the pillowcase is not recommended.

3/ Keep the device in its packaging in a dry place, away from direct sunlight.

4/ Protect the device from high temperatures, open flames, extreme temperatures, humidity and impact.

5/ Store the product at temperatures between 0-40°C.