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
INSTRUCTION MANUAL MASSAGE PILLOW CHOICE Dr. **DC-022**

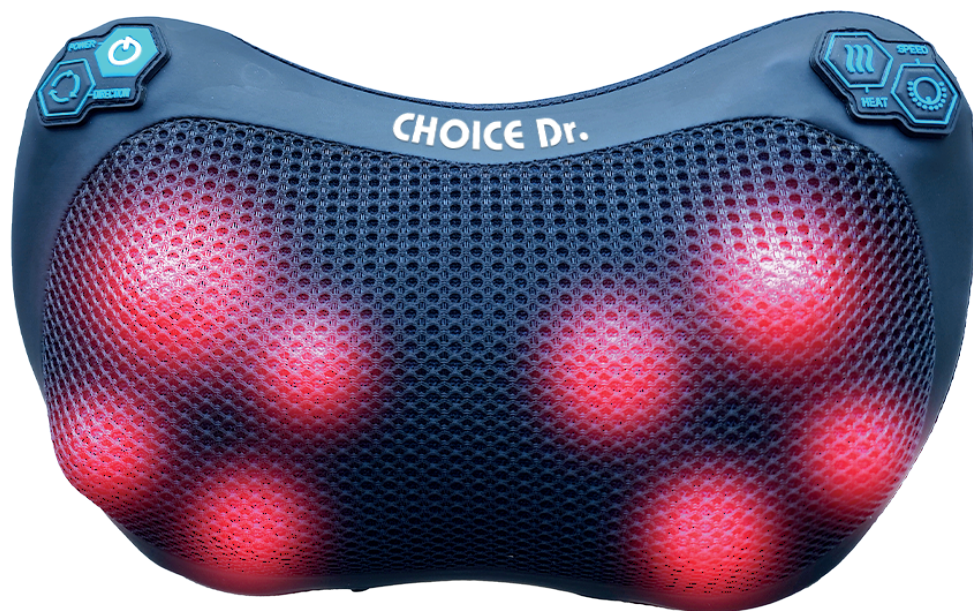CHOICE Dr.

Thank you for choosing and supporting CHOICE Dr.

We hope that you are satisfied with our product and service, and that we have the opportunity to serve and accompany you and your family.

If you have any feedback or reviews regarding the quality of our product and service, please contact our hotline at **1800 8379** so that we can receive and promptly address them.

In addition, please scan the **Zalo OA** code to activate the electronic warranty, receive additional gifts, and enjoy promotion and **24/7** support from CHOICE Dr.

Thank you once again and we are happy to have take care of you and your family's health!

Choice Dr

CONTENT

NOTE

HOW TO USE THE MASSAGE PILLOW

PRODUCT FEATURES

NOTE

PLEASE DO NOT USE THE MASSAGE PILLOW OR ANY VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY:

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device but we advise you to consult a doctor first.

Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in a number of cases. This must be done on the advice of and in accompaniment with a doctor, specialist or physiotherapist.

HOW TO USE THE MASSAGE PILLOW

USER MANUAL

1. To use the massage pillow, please connect to **home charger (6)** hoặc **or car charger (7)** come with the package.
2. After connect to charger, short press **button (1)** to start.
3. To **change** massage head **rotation direction** , please press **button (2)**.
4. To **change massage speed**, please press **button (4)**, there are total 3 level massage speed.
5. By press **button (3)** you can turn on/off **heat function**.

PRODUCT FEATURES

- 1 Three-stage speed regulation
- 2 Heating mode - infratherapy
- 3 Changing the direction of the massage
- 4 Balls that massage pillows located under a special fabric that allows the skin to breathe during the massage
- 5 Built-in timer that turns off the device after 15 minutes of operation
- 6 Car and traditional power supply
- 7 Intuitive buttons to control the massager
- 8 Headrest elastic band

SPECIFICATIONS

Product code: DC-022	Rated working time: 15 minutes
Rated Voltage: AC220V/50HZ	Net weight: 2,1kg
Rated Power: 40W	Color box size: 33*11*19 cm