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Model: **MAGIC ECO PRO-81**

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PORTABLE WATER FLOSSER

Model: PRO-81

Thank you for choosing and supporting MAGIC Vietnam!


We hope that you are satisfied with our product and service, and that we have the opportunity to serve and accompany you and our family.

If you have any feedback or reviews regarding the quality of our product and service, please contact our hotline at **1800 8379** so that we can receive and promptly address them.

In addition, please scan the **Zalo OA** code to activate the electronic warranty, receive additional gifts, and enjoy **6 months of warranty** and **24/7** support from MAGIC.

Once again, Magic would like to thank you and we are happy to help you take care of your oral health and hygiene.

PRODUCT INTRODUCTION

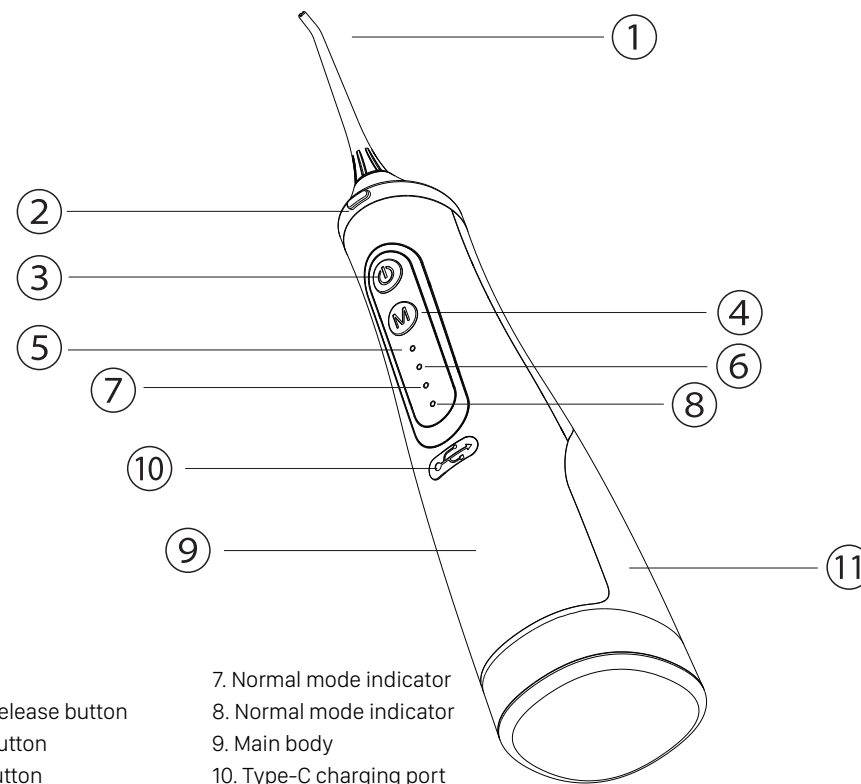- | | |
|--------------------------|--------------------------|
| 1. Nozzle | 7. Normal mode indicator |
| 2. Nozzle release button | 8. Normal mode indicator |
| 3. Power button | 9. Main body |
| 4. Mode button | 10. Type-C charging port |
| 5. Pulse mode indicator | 11. Water tank |
| 6. Gentle mode indicator | |

PRODUCT BASIC PARAMETERS


- | | |
|---|--|
| 1. Water Tank: 310ml | 8. Noise value: <60 db |
| 2. Water Pressure: 30-120PSI | 9. Certificates: CE, RoHS |
| 3. Battery: 1500mAh 5C power battery | 10. Memory mode: YES |
| 4. Waterproof Level : IPX7 | 11. 360° rotating Nozzle |
| 5. Material: ABS&PC | 12. Product size: 241*84*70 mm |
| 6. Power Supply: Waterproof Type-C | 13. Accessories: 1 machine, 4 nozzles, |
| 7. Working Model : Pulse, Soft, Normal, Power | 1 charging cable, 1 manual |

PRODUCT MAINTENANCE


1. Please use water or neutral detergent to clean the product. Do not use corrosive or abrasive detergents as these can damage the product. Do not use hot water above 40°C.
2. Do not immerse the main body into water, wipe clean with a dry or wet cloth, or clean it with running water and dry it with a rag.
3. It is recommended that the nozzles are replaced every 3 months.

INSTALLATION AND USAGE



1. Insert the nozzle into the nozzle socket, hear the "click" sound to indicate the installation in right place


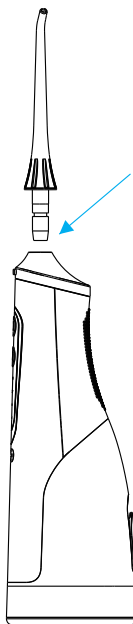
2. Open the water tank cover and fill the tank with water, do not fill with warm water above 40°C, salt, ointment or chemicals. After filling with water, please close the water tank cover

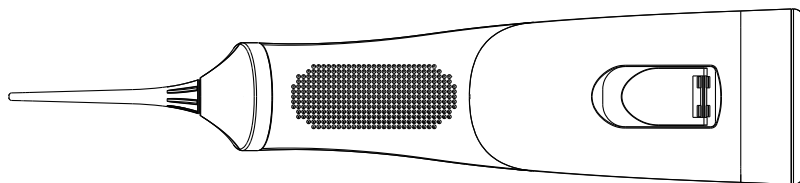3. Tap the mode button, select the mode, put the nozzle into the mouth and press the power button to turn on the machine for cleaning teeth operation.

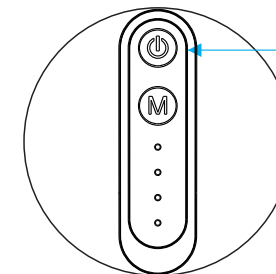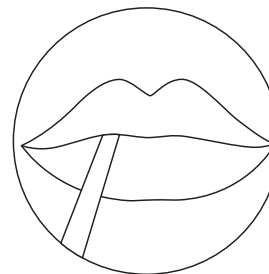4. When using, lean forward slightly, open your mouth slightly and let the water flow out naturally from the corner of your mouth.
Do not stand or sit to flush your teeth, and then do not open your mouth widely to flush your teeth to avoid splashing water on your body or face.
For best results, start from the molars (back teeth) and work towards the incisors.
When cleaning along the gum line, pause briefly between the teeth to ensure that the area around the gums and between the teeth is cleaned.

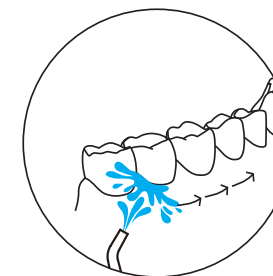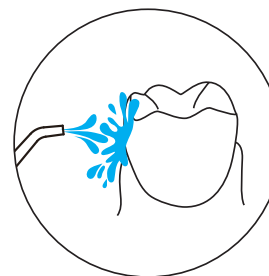5. Mode selection


The water dental flosser is switched to the pulse mode in default state.

The flosser cycles through the four modes: pulse, gentle, normal and strong in a top-to-bottom sequence. It is recommended to use the gentle mode for the first time.


6. Charging guide

Open the Type-C charging cable and plug one end firmly into the charging port of the product and the other end firmly into the power adapter to charge.


If there is water inside the charging port of the machine, at first wipe it off before plugging in the type-C charging cable.

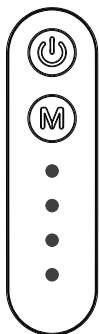The four indicator light flash at the same time, signifying that the battery power is low.

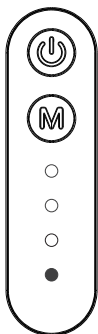The four indicator light flash one by one, signifying it is being charged.

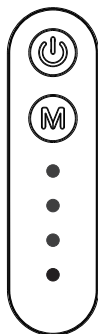The four indicator light without any flash, signifying that the battery is fully charged.

WARNINGS

Use and Storage

1. Don't dismantle the product or change the parts of the product, so as not to damage the product or cause electric shock. If you did, the warranty is invalid.
2. Don't use water above 40°C to fill the container to clean the mouth to avoid burns.
3. If the charging product into the water, at first unplug the power adapter on the power socket, and then fish up the product to avoid the risk of electric shock.
4. Don't put the product into the fire or heating, do not charge it at high temperatures or humidity.
5. Do not put the product in a humid place for a long time, and avoid the direct sunlight.
6. Don't use wet hands to pull out the power plug and charging plug during the charging process.
7. This product does not come with a power adapter, the output voltage of customer's own adapter must choose DC 5V, such as cell phone power adapter, power adapter before use, please check whether the voltage meets the voltage requirements of the local use.
8. If the product is not used in a long time, please unplug the power cord and charging plug to avoid fire or other dangers.

COMMON FAULTS AND TREATMENTS

PHENOMENON	POSSIBLE CAUSES	SOLUTION
Water tank leakage	Water tank has a small air inlet hole	The air inlet hole is to keep the air pressure inside the water tank in balance with the outside air pressure, the machine can be used vertically
The water pressure is too low	The power is not enough	Full charge and then use
	The air inlet hole of the water tank is blocked	Clean the air inlet hole or open the water tank filler cap
After charging, the machine can only work for a few minutes	Battery life is exhausted (about 500 times)	Contact your dealer
	Faulty charging cable or charging head	Replace charging head or cable
There is a sound without water	Internal valve cocked	Add water through the nozzle filler hole under the power on