

POONGSAN
Premium

Luxury Massage Chair

USER MANUAL

**PLEASE KEEP THIS MANUAL
FOR FURTHER USE**



MCP-902


Please read and follow all safety and operating instructions listed on this user manual before use. In order to use it in a right way, please keep the user manual properly for future reference.

Safety Precautions	2
Product Specifications	4
Product Components	5
Chair Installation	7
Before Use	9
Operation Guide	10
Maintenance	17
Trouble Shooting	18


If the sounds or movements are not normal, please turn off the massage chair and remove the plug from the electric outlet. Please contact the customer service department for further assistance.

Problem	Solution
The massage chair is connected to power but is not functioning.	<ol style="list-style-type: none"> 1. Check to see if the power plug is inserted properly or if the power outlet is functioning properly. 2. Check to see if the main power switch is in the "I" position.
The remote controller is responding, but the massage chair is not functioning.	The massage chair may have automatically shut down due to overheating protection devices. Please turn off the power and allow the chair to cool off for 30-50 minutes.
Both the remote controller and massage chair are not responding.	The massage chair automatically turns into "full stop" or into "idle" mode when a massage program is completed. Please press the power button again to turn on the chair.

MAINTENANCE

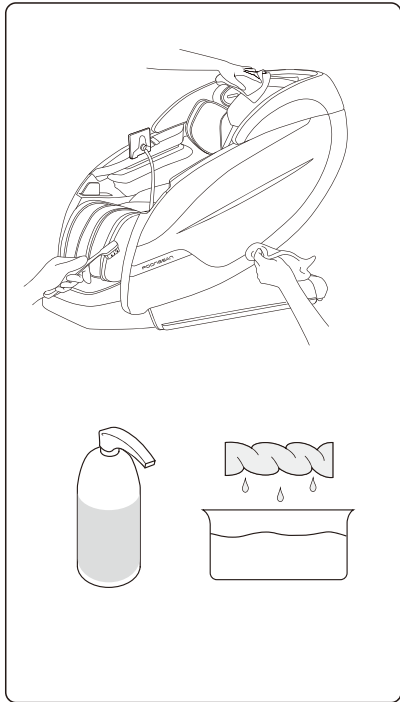

- Unplug the power plug before cleaning.
- Do not pour any liquids directly to the chair. Spray or dampen the cleaning cloth.
- Use a non- toxic cleaner and dry cloth to clean the chair.
- Be sure the chair is completely dry before plugging the power plug into the outlet.
- Do not use any drying machine to dry the areas cleaned. Only air dry to avoid any damage to the chair.
- Do not use any product which contain benzene or alcohol to clean the chair.

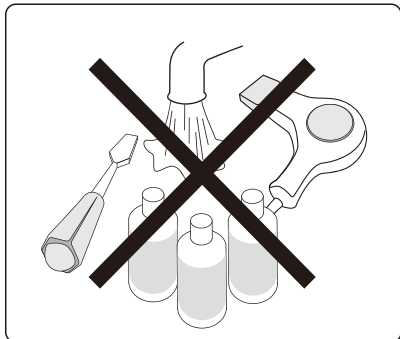

SAFETY PRECAUTIONS

Warning!

Prior to use the massage chair please read all the operational and safety guidelines.

The power plug must be fully inserted into an appropriate outlet. Do not modify the plug or remove the ground wire which may damage the massage chair.

1. Make sure that the power plug is securely inserted into the proper outlet to reduce the risk of short circuits and fires.
2. Pull out the power plug from the massage chair before cleaning or after using.
3. Please do not use when the humidity is high.
4. Do not place the massage chair near or in areas with high temperature.
5. Do not use the massage chair if the power cord is damaged.
6. For those people who suffering from illness, please consult doctor before use.
7. The massage chair is not suitable for children aged 12 and below.
8. Do not sleep on the massage chair while using to avoid any physical injury.
9. Do not use the massage chair more than 30 minutes a day. Monitor your condition daily for any muscle or body pain or discomfort please refrain from using the chair and see your doctor.
10. Do not move the massage chair by pulling the power supply cord.
11. Do not stand on the massage chair and its legrest.
12. When moving the massage chair please follow the instructions to avoid structural damage.
13. When using the massage chair please make sure that there are no children or pets nearby, behind or front when the massage chair is in motion it may cause serious injury.
14. Do not modify your massage chair which will void your warranty. Please contact the authorize technician to assist you.
15. This massage chair is for personal use and designed for household use only.

SAFETY PRECAUTIONS

16. Always check whether the surface of the massage chair is torn or ripped before use, if you find any unfamiliar rips in any area of the cloth of the backrest, do not use the massage chair. Kindly contact the service center immediately for repair.
17. Please do not use if you have a heart condition or have a pacemaker.
18. If you experience pain during usage, stop using the massage chair immediately, if any injury occur seek medical treatment immediately.
19. Do not massage the swollen or inflamed parts of the body.
20. Do not use this product as substitute for medical treatment.
21. Each massage session should not exceed 30 minutes, massage the targeted area only for 15 minutes, over massaging can lead to overheating, thus giving the appliance a shorter lifespan.
22. Do not operate where the aerosol (spray) products are used or when oxygen is being administered.

NOTE:

This appliance is not intended for use by person (including children) to reduce the physical sensory or mental capabilities or lack of knowledge how to use the massage chair.

Children should be supervised to ensure that they do not play with the massage chair.

If the power supply cord is damaged, it must be replaced by a special or assembly available from the manufacture or its service agent.

The massage chair heats up during operation, persons sensitive to heat must take extra caution when using the massage chair.

OPERATION GUIDE



The method of playing the bluetooth music in massage chair :

A: Match: Open the bluetooth function in mobile phone or Ipad ,and search nearby. Choose the massage chair bluetooth music equipment and match with it.

(Note: The massage chair Blue tooth code MCP-902)

B: Play music:After matching with it, you can display the music by mobile or APP music display software.


After use

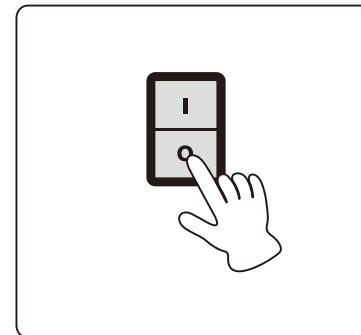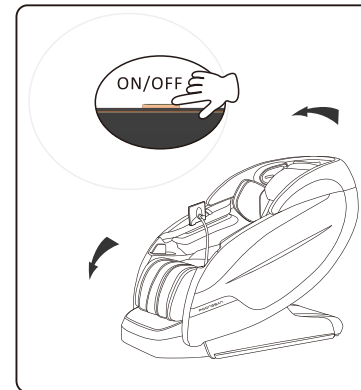

1. After the massage program is completed, please press the switch button, the massage chair will be automatically recover.
 - Please check surrounding area prior to adjust the recline angle, pay special attention to children or pets that may be close proximity.
 - Please make sure your arms and legs are away from gaps between the armrest, legrest, and footrest, to prevent any possible pinching or injury.
2. Press the power button to power off the massage chair and return remote controller to the storage box.
3. Turn off the main power switch.
4. Remove plug from power outlet.
 - * Please do not pull on the power cord to remove the plug.
 - * Please do not unplug with wet hands.

OPERATION GUIDE


7. Adjust seat position

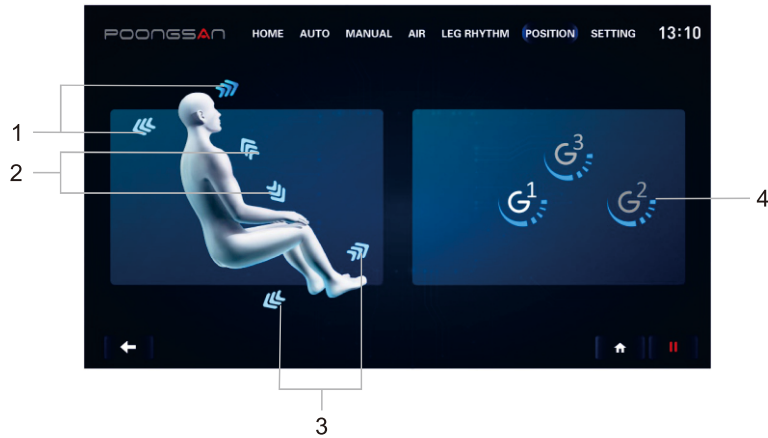

1. Incline and raise the backrest
2. Stretch legs and thighs
3. Lift and lower thighs and legs
4. Zero gravity: 3 boxes to choose

Note: Press relevant keys on the touch screen for adjustment.

8. Set up

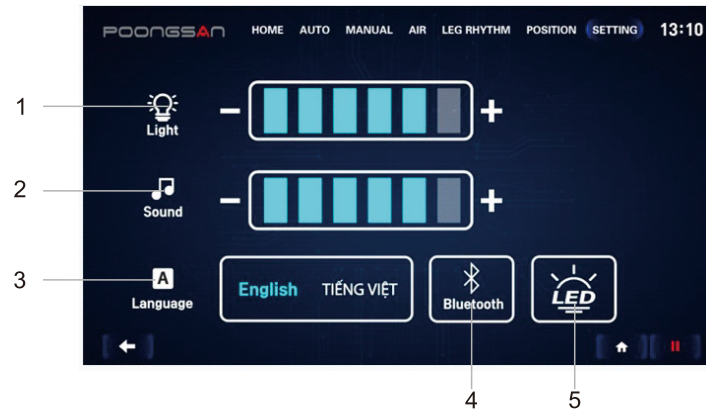


1. Adjust brightness
2. Adjust sound volume, there are 6 levels volume to adjust bluetooth music and voice controller.
3. Choose language
4. Bluetooth, note: make sure bluetooth is open when use the voice control function.
5. LED light: Press to the symbol to on/off the LED light on the speaker position.


PRODUCT SPECIFICATIONS

Description	Specifications
Product Name	Massage Chair
Model No.	MCP-902
Rated Voltage	220-240V~
Rated Frequency	50/60Hz
Rated Power Consumption	230W Max.
Dimension (LxWxH)	Main Body: 1460x856x1245mm
	Legrest: 600x535x630mm
Weight	Net Weight: 135kg+25kg
	Gross Weight: 156kg+28kg

Main Body Diagram (Front)

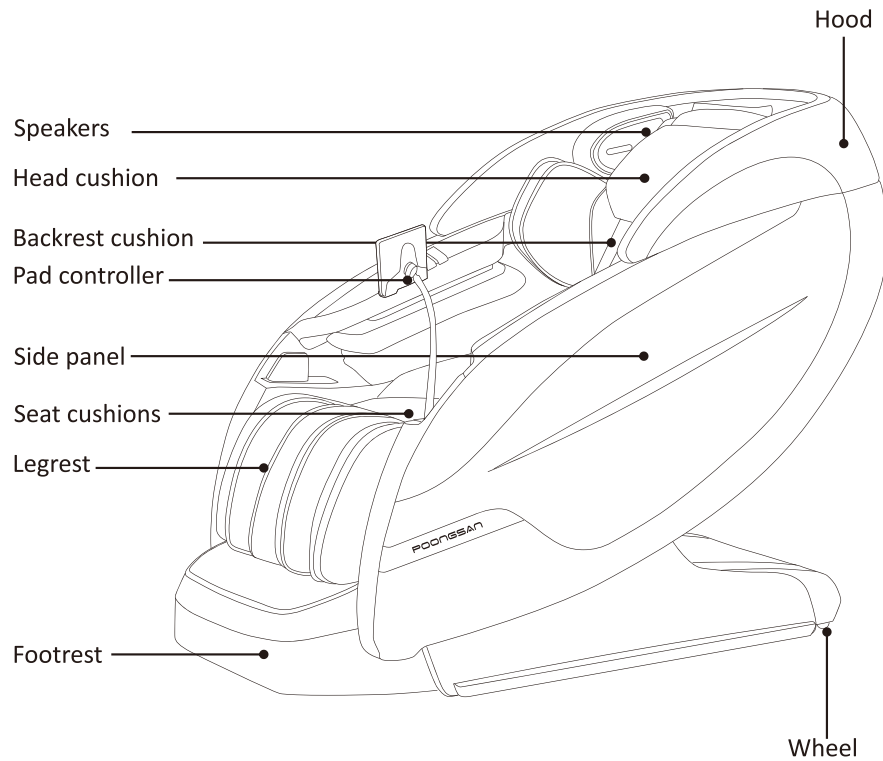


5. Airbag

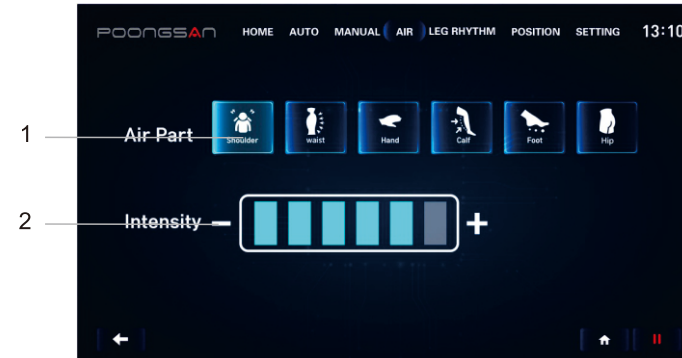

1. Positions of airbags(6 positions to choose)
 2. Adjust airbag intensity.
- Note:** Press relevant keys on the touch screen for adjustment.

6. Leg vibration

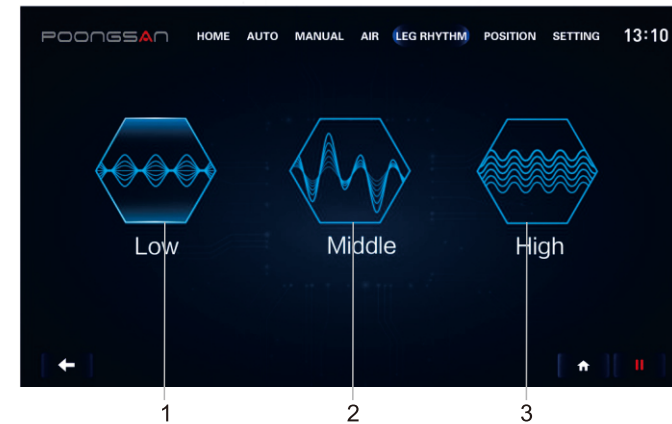

- 1.Slow speed vibration
- 2.Medium speed vibration
- 3.High speed vibration

OPERATION GUIDE



- 1. 4D massage intensity of massage chair
- 2. Massage intensity of airbag
- 3. Massage power of massage chair
- 4. Massage time
- 5. Enter auto massage menu to select interface
- 6. Warming upper half of body
- 7. Warming feet
- 8. Mechanism roller heating
- 9. Rolling and pressing soles
- 10. Rolling and pressing thighs
- 11. Rolling and pressing arms.

Note: Press relevant keys on the touch screen for adjustment.



4. Manual mode



- 1. Shiatsu
- 2. Tapping
- 3. Kneading
- 4. Rolling
- 5. Swish
- 6. Knocking
- 7. Select up/down
- 8. Select massage point
- 9. Adjust massage zone (wide, medium, narrow)
- 10. Select massage zone
- 11. Full massage

PRODUCT COMPONENTS

Main Body Diagram (Rear)

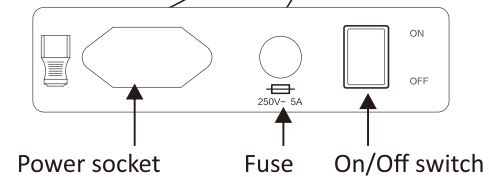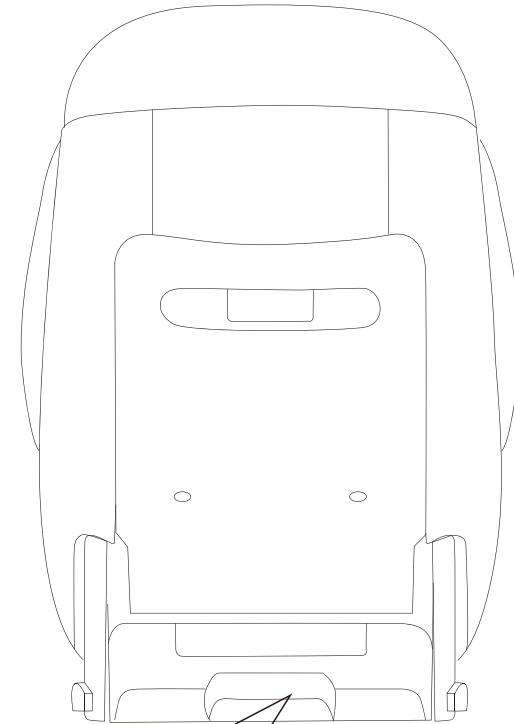


Power socket Fuse On/Off switch

Installation Method


Step 1: Installation of legrest unit

- ▶ Firstly connect the air-pipes and wires of legrest with the mainbody of massage chair.
- ▶ Lift the leg and insert the round hole in the upper right corner of the leg into the metal pin shaft of the massage chair, and then insert it with a buckle to fix it.
- ▶ When installing the left side, first move the whole leg to the right, align the round hole in the upper left corner and the pin of the massage chair, move the lower leg to the left, and then insert the buckle to fix it.
- ▶ Pull the holster on the leg to the zipper of the massage chair.
- ▶ Finish the leather case.

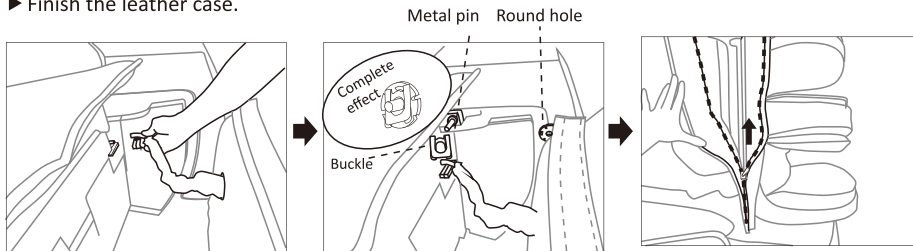


Step 2: Installation of controller unit.

- ▶ Put the pad controller bracket into the left slot of the massage chair first.
- ▶ Then dock the pad controller cable with the connector on the massage chair.

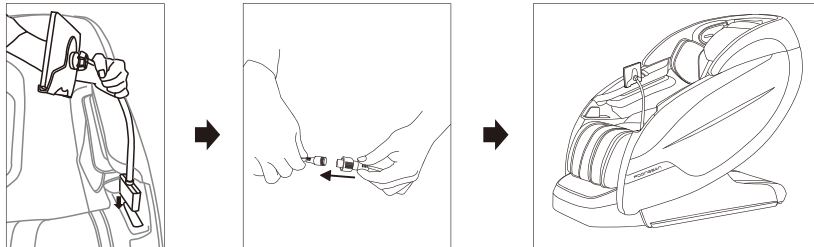

2. Manual control panel interface and description

-Open manual control mode, control menu shall display. Choose menu that you want.

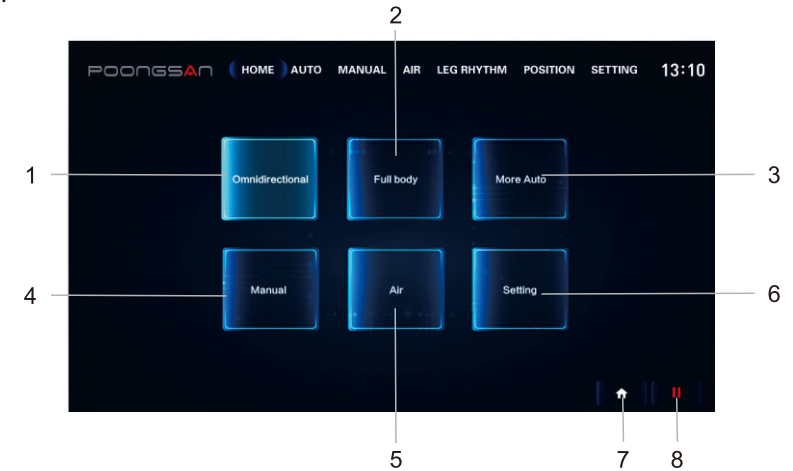

- 1.Auto mode (Omnidirectional) 2.Auto mode(Full body massage)
 3.More auto 4.Manual 5.Airbag massage 6.Set-up
 7.Main menu 8.Pause



3. Auto mode



There are 12 auto programs in total.

Note: The default time for each auto program is 15 minutes.


1.Body scanning

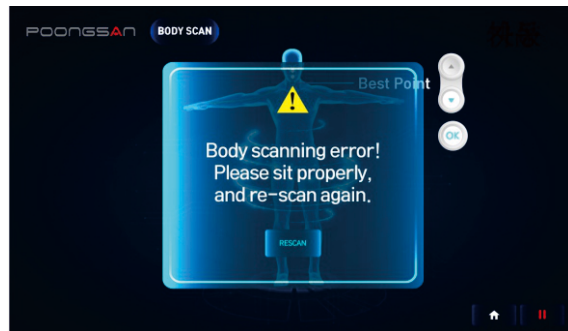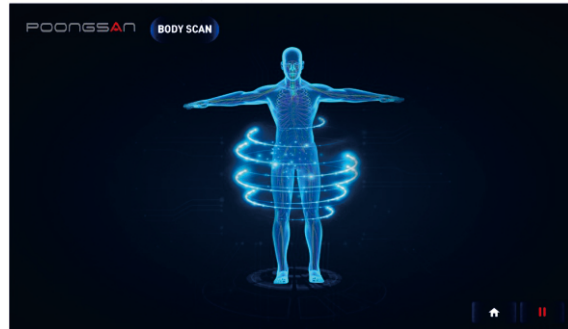

- 1.Latest intelligent body position detection.
- 2.After the detection is completed, the user can adjust the shoulder position according to the selection within 8 seconds or press OK directly to confirm the shoulder position.
- 3.When the shoulder detection fails, the prompt page will pop up. At this time, please adjust your sitting position and press "RESCAN" to conduct human detection again.

How to move


As shown in illustration, please do not use the armrest when moving the massage chair, but hold the bottom of the product and the backrest for movement.

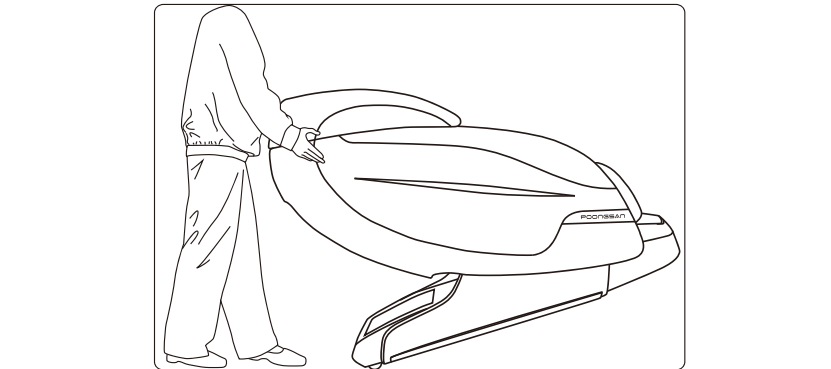

Caution:


- 1.Do not move the massage chair when it is occupied.
- 2.Unplug the massage chair cable from power source.
- 3.Do not drop or suddenly let go of the chair when moving as it may damage the chair.
- 4.Two person participation is recommended.
- 5.Placing a mat underneath the massage chair is advisable, as moving process might induce scratches on the floor.

BEFORE USE

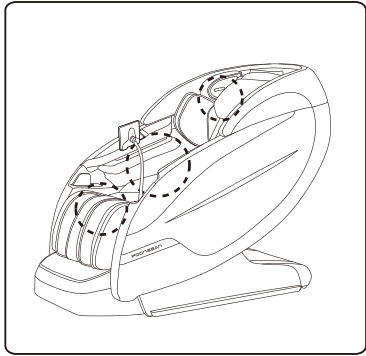


1. Make sure that the product is in good condition

- Make sure that headrest and backrest pad are normal.
- Make sure that the power supply cord is normal.
- Make sure whether the footrest is erect. (if it is under the horizontal surface, when you sit on the massage chair and exert big force the chair will be damage).

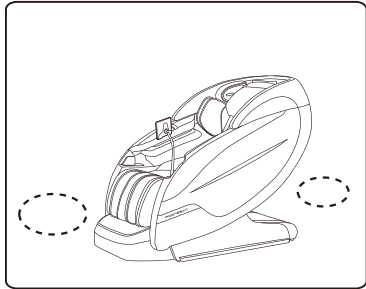


2. Check the surrounding environment

- Make sure there is no barriers such as baby and pets around the massage chair.

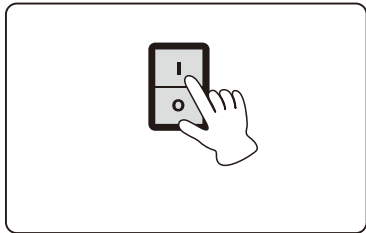

3. Connect the power supply cord of the massage chair

- Connect the power supply cord.
Turn OFF the power switch on the back of the massage chair

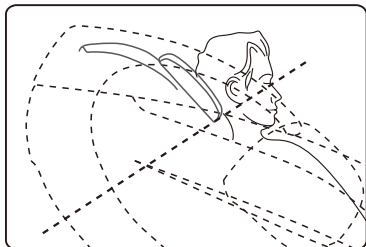

4. Adjusting the position of the headrest when using

- Adjust the under-part position of the headrest and make the position of the ear and on the same horizontal plane.

OPERATION GUIDE

Pad Controller

