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Model: **MAGIC ECO PRO-82**

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PORTABLE WATER FLOSSER

Model: PRO-82

Thank you for choosing and supporting MAGIC Vietnam!

We hope that you are satisfied with our product and service, and that we have the opportunity to serve and accompany you and your family.

If you have any feedback or reviews regarding the quality of our product and service, please contact our hotline at **1800 8379** so that we can receive and promptly address them.

In addition, please scan the **Zalo OA** code to activate the electronic warranty, receive additional gifts, and enjoy promotion and **24/7** support from MAGIC.

Once again, Magic would like to thank you and we are happy to help you take care of your oral health and hygiene.

SPECIFICATIONS


1. Container: 350ml
2. Voltage: 5V~2A
3. Power: 3W
4. Battery: 1800mAh
5. Charging time: 4 hours
6. Pressure: 30 - 150 psi
7. Noise level: less than 72dB
8. Weight: 334g
9. Dimensions: 80 x 80 x 312 mm

WARNING


To reduce the risk of injury to persons:

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Oral irrigators shall be used only with cold or warm water, or specific solutions as defined in the manufacturer's instructions.
- Children shall not play with the appliance. It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- The product contains batteries that are non-replaceable.
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charge.
- If product is opened/disassembled for ANY reason, warranty is VOID.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Please contact your local distributor.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba *Naegleria fowleri* may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet.
- Do not drop or insert any foreign object into any opening on the product.
- Keep charging unit and cord away from heated surfaces.
- Do not use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Children should be supervised to ensure that they do not play with the appliance.


PRODUCT DESCRIPTION

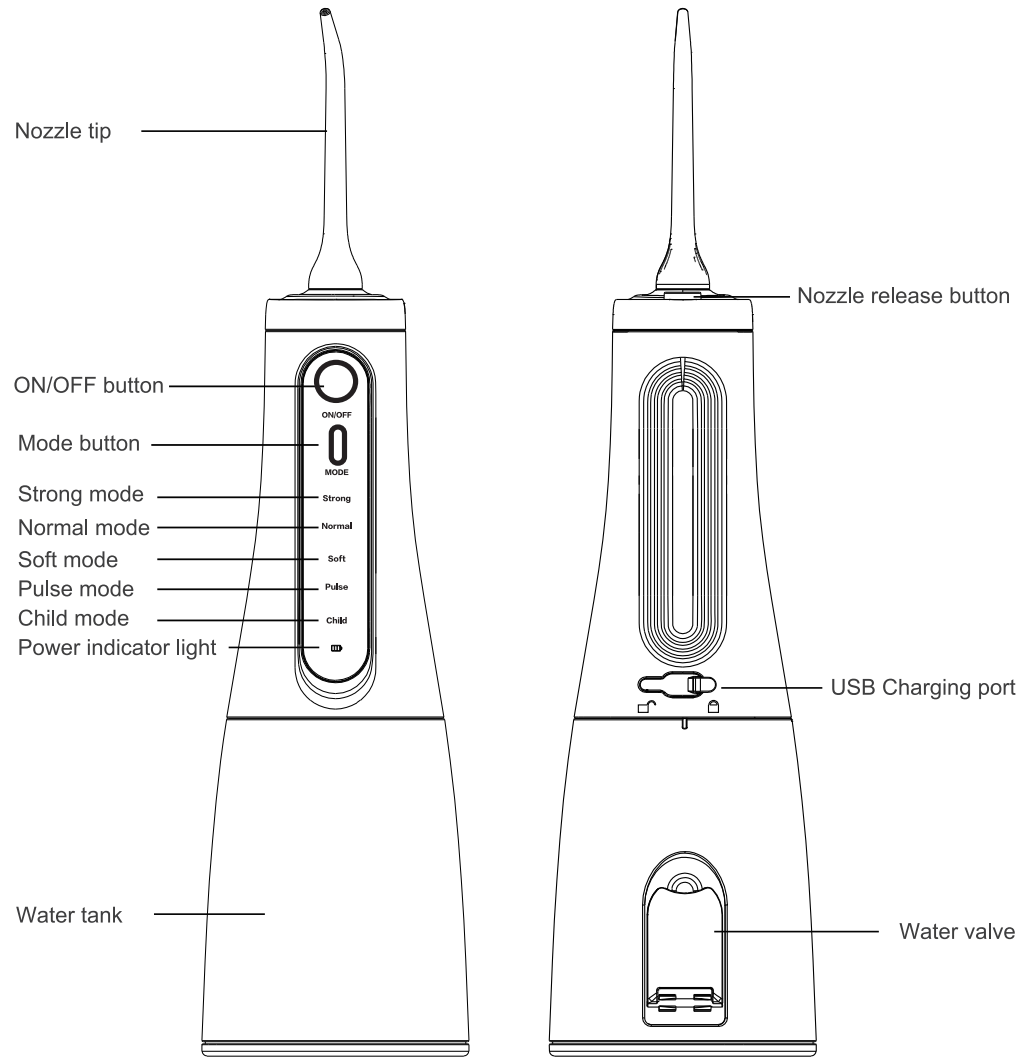PRODUCT ACCESSORIES

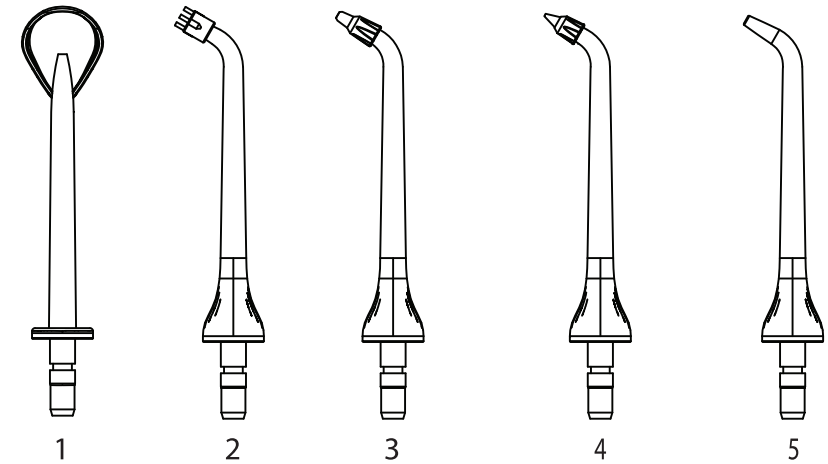1. Tongue cleaner *t1
2. Dental plaque tip *1
3. Orthodontics tip *1
4. Periodontal pocket tip *1
5. Classic jet tip *2

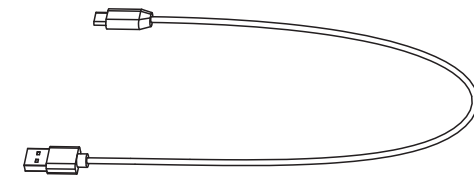5. USB charging cable *1

The jet tips and other changeable components can be purchased from retail stores

OPERATING INSTRUCTION


1. Preparation before starting

Very important! For the first time to use, please fully charge the device. During charging, the product CANNOT be used.


3. Pour water into the water tank

Note: Please ensure that when pouring water into the water tank, the device is switched off.

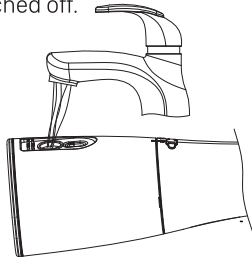
5. Power on

Put the nozzle into your mouth and press the ON/OFF button to start using.

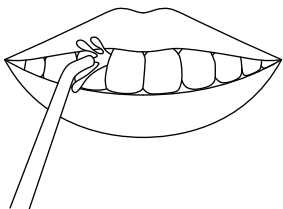2. Insert nozzle

Insert the nozzle into the hole until you hear the "click" sound.

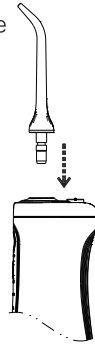
4. Mode selection

Press the MODE button to select your suitable mode.

6. Change nozzle



Power off and remove the tip, press the nozzle release button, and then pull up the nozzle.

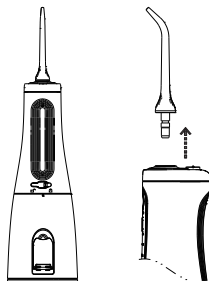- ON/OFF Button is the key to power on and off. The first short press is to start up, the second short press is to power off.
- Mode Button is the working mode selection button. Short press to switch the 5 working modes: STRONG-NORMAL-SOFT-PULSE-CHILD. We can select the mode firstly before powering on the device, or select the mode while using.
- The device CANNOT be used during charging.
- The charging indicator light keeps red on when being charging, and keeps white light on when fully charged. If the charging indicator light becomes red and flashes, it means the battery is required to recharge.

USAGE TIP

- Tilt down to place the nozzle in the mouth, and point the nozzle at the tooth, then turn on the device.

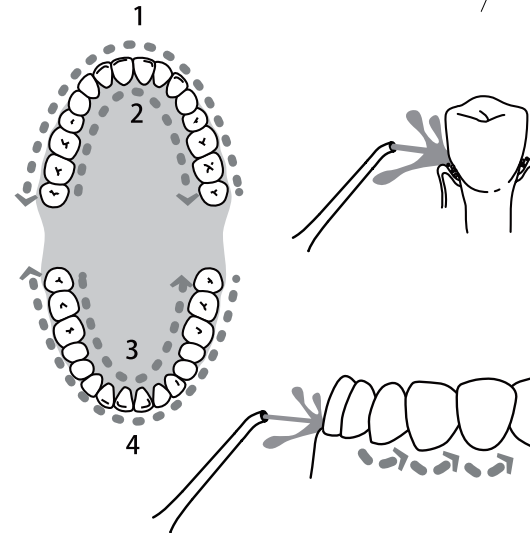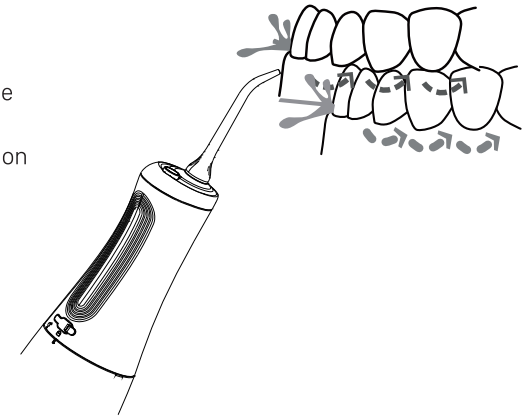
- Make water jet at 90° to the gum, and close your lips slightly to avoid water spray. For best results, start with the molars (back-teeth) and clean towards the front teeth. Slide the tip along the gum line and keep flossing a while at interdental places, till everywhere are cleaned.

ASSEMBLY AND DISASSEMBLY


Clean water tank

Turn off the device, and open the back tank lid with the injection hole facing down, then turn on the device after clearing the residual liquid in the tank, to drain the residual liquid in the device until there is no more water to drain, then turn off the device.

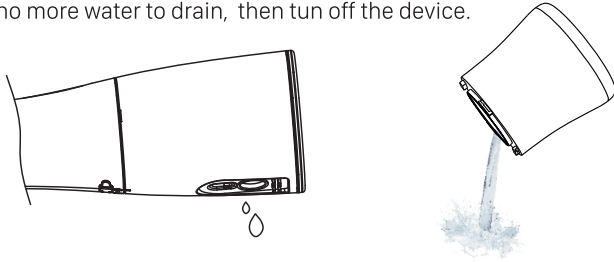Remove and clean

- Using both hands as indicated in Picture below
- Keep your left hand to hold the product body, and your right hand to hold the water tank and take off it.
- You can use brushes to clean water tank if necessary. (Note: brushes are not included).

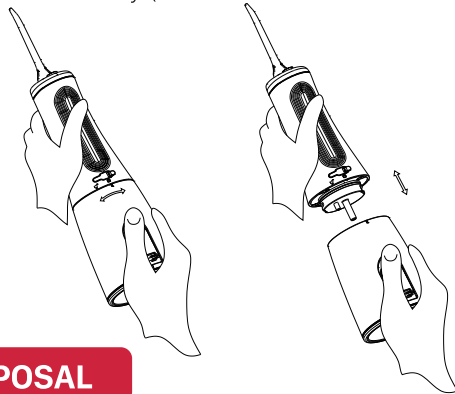BATTERY CARE AND DISPOSAL

Battery Maintenance

- For the first use, the device must be fully charged approximately 5 hours.
- If the device is stored for a long time, please fully charge it in advance.
- In order to extend battery lifetime, please recharge it in time when power is insufficient.

Battery Care and Disposal

- The device includes a lithium battery which cannot be changed. To preserve the environment, do not throw away the appliance with the normal household waste at the end of its life, hand it over to local recycling center

FREQUENTLY ASKED QUESTIONS

Problems	Reasons	Solutions
Leakage of water between jet tip and screw part	The seal ring on the fixed socket for the nozzle is damaged	Contact the after-sales
	The nozzle does not install correctly	Re-install correctly
Water tank leaks water	Seal ring is damaged	Replace the seal ring, call after-sales
Inadequate pressure	Nozzles damaged	Replace with new nozzle
	Low voltage	Charging
	Low battery	Charging
Water is not coming out	Nozzle blocking	Suck the jet tip by mouth
	Water pipe blocking	Inject the water from the nozzle
		Shake after fill the water
Short working time after full charging	Insufficient charging time	Charging
	Circuit fault	Contact after-sales

If the above methods cannot resolve these problems, please consult with the retailer or proceed to the designated after-sales repair site for repairs.

Gentle Reminder

1. For the first-time usage, please select "Child mode", in order to avoid the gums bleeding due to the lack of adaption to usage.
2. This product is not suitable for children, infants and other operators to use, unless they get their safety responsible person in the use of this product supervision or guidance. Keep out of reach of children or infants. Do not use as toys
3. For patients with periodontal disease or those who have had oral surgery in the last 2 months, please consult your doctor before use.