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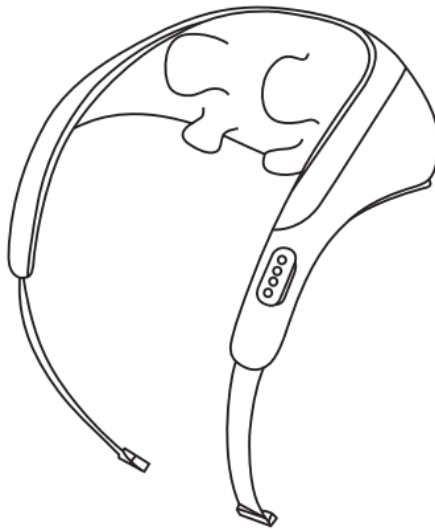
INSTRUCTION MANUAL

NECK AND SHOULDER MASSAGER


MODEL: MK395

Read this instruction manual carefully before using

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
Product Instructions

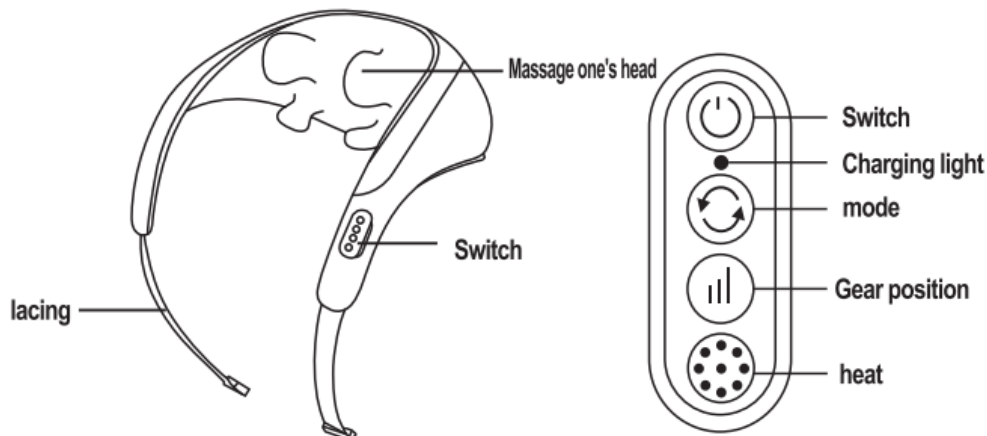Product Parameters

Product name	Rated voltage	Rated power	Product size	Product weight
Neck and Shoulder Massager	DC 5V	25 W	≈ 23 x 13.5 x 16 cm	≈ 1.35 Kg

Correct usage method

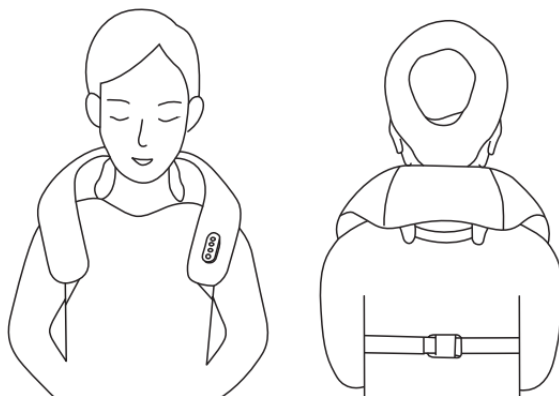*** Preparation before use**

Please check if the product and accessories are complete and undamaged.

Press and hold the “switch” button for 2 seconds, and the product starts working.

***Mode selection:**

The default mode for startup is mode 1. Short press the “mode” button to switch between mode 1 and 2.

***Heating mode**

Boot the default heating off, short press the “heat” button, heating high and heating low in turn and turn off the heating, so as to cycle.

***Strength mode**

The default mode is second gear when starting up. Short press the gear button to cycle the first to the third gear to second gear.

***Product shutdown**

Automatic shutdown: After a total of 10 minutes of use, the product stops working.

Manual shutdown: Press the “switch” button to stop working.


Charging method

1. This product cannot be used while charging.
2. Plug the charging cable into the charging port on the product.
3. When the product is working, it is connected to the power supply to charge, and the product starts the self-protection program Dynamic stop working.
4. The product charging indicator light is red, and after charging is completed, the indicator light is green.

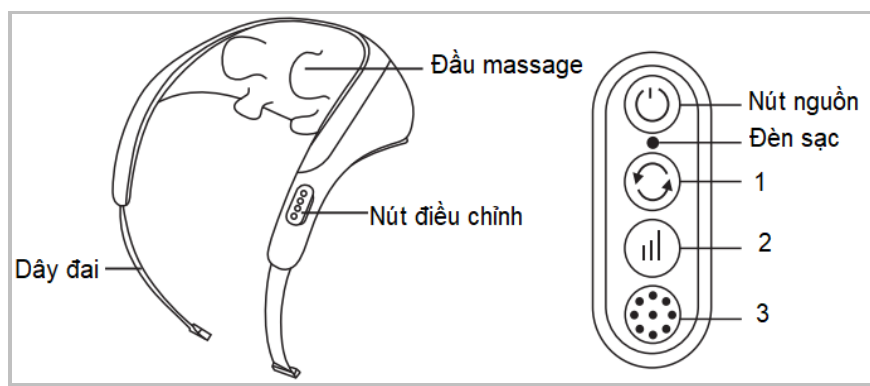* It is normal for the main and power adapter to heat up slightly during charging.

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