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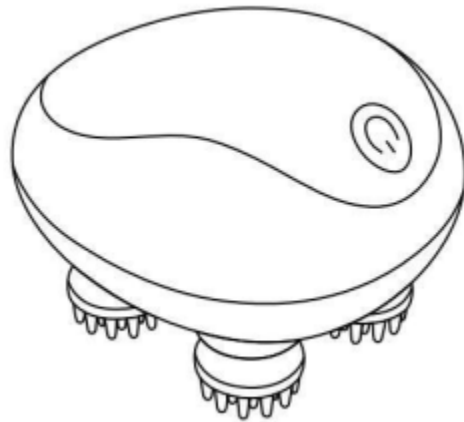
# INSTRUCTION MANUAL

## HEAD MASSAGER

MODEL: MK399

Read this instruction manual carefully before using

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**ENG**

## **Product Parameters**

Product name	Rated voltage	Rated power	Product size	Product weight
Head Massager, Model MK399	DC 5V	5 W	11.5 x 7.3 x 9.5 cm	192 Gram

### **Note:**

1. Do not use it without authorization if the head has been operated on.
2. Do not use it if the head has skin damage, swelling, inflammation, congestion and other symptoms.
3. Do not use the massager without permission for people with slow motion, physical and mental maladjustment or intellectual impairment.
4. The massager is designed as a household electronic health care product, which is not suitable for diagnosis, treatment and other commercial purposes.
5. The massager is a personal care product, which can not replace the advice and treatment provided by medical professionals.
6. Children and pregnant women are not allowed to use the massager.
7. Do not use the massager while driving. Follow the safe driving rules.
8. Do not hit, fall or use sharp objects to scratch the massager.
9. Minors must use the massager under the guidance of their parents or persons with supervision ability.
10. Stop using the massager immediately if you start to use it or if you are not in good time.
11. Do not use the massager after drinking or when you are not in good time

### **Correct usage method**

When using the product, hold the product key with one hand, put the massage head of the host on the head and massage each position of the head

1. Long press the on key of the main machine for two seconds, the main machine starts and the massager starts to work in three modes

Mode 1: massage claw forward slow rotation

Mode 2: press slow rotation alternately

Mode 3: the massage claw rotates alternately in the positive direction

2. Indicator light Description: red and green light, the product can not work when charging

Green light--The green light is always on when the product is working;

Red light--Flashing when charging and low power;

When charging--The red light is flashing;

After full charge--The green light is always on;

### **Charging method**

1. The built-in lithium battery of this product can be recharged repeatedly. Please charge in time after weak current automatic shutdown protection
2. In the process of using the massager, if the indicator light flashes red, it will remind that the power is low, the massage speed or intensity is relatively weak (not machine fault), and it should be charged in time
3. Please use 5vm1a charger to charge the machine. Do not use the charger when charging
4. If the rechargeable battery is not used for a long time, the service life of the battery will be affected. It is recommended to charge the battery at least once every 3 months.

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